easyGWAS : An integrated interspecies platform for performing genome-wide association studies

Dominik Grimm¹, Bastian Greshake¹, Stefan Kleeberger¹, Christoph Lippert¹, Oliver Stegle¹, Bernhard Schölkopf², Detlef Weigel³ and Karsten Borgwardt^{1,4}

¹Machine Learning and Computational Biology Research Group, Max Planck Institute for Developmental Biology and Max Planck Institute for Intelligent Systems, Tübingen, Germany

²Department of Empirical Inference, Max Planck Institute for Intelligent Systems, Tübingen, Germany

³Department of Molecular Biology, Max Planck Institute for Developmental Biology, Tübingen, Germany

⁴Center for Bioinformatics Tübingen, Eberhard Karls Universität, Tübingen, Germany

Contents

1	1 Motivation		3
	1.1 Difficulties in performing GWAS across traits and species .		3
	1.2 Role of $easyGWAS$		4
	1.3 Functionality \ldots \ldots \ldots \ldots \ldots \ldots \ldots \ldots		4
	1.4 Inter-compatibility with statistical genetics software package	25	4
2	2 Example		4
3	3 Conclusion and future plans		6
\mathbf{A}	A Data: Genotypic & phenotypic data and meta informatio	on	7
	A.1 Available published data		7
	A.2 How to upload new data \ldots \ldots \ldots \ldots \ldots \ldots \ldots		8
	A.2.1 Phenotypic data and meta information		8
	A.2.2 New genotypic data		8
в	B Integrated methods to perform genome-wide association	studies	8
	B.1 Methods to perform a GWAS		9
	B.2 Transformations to standardize data		10
С	C The web application interface		10
	C.1 The $easyGWAS$ wizard and experiment history		10
	C.2. The data center		13

	C.3	The download center	14
D	The	web application backend	14
\mathbf{E}	Tut	orials	15
	E.1	How to perform a GWAS easily?	15
	E.2	How to upload an own phenotype and perform a GWAS on it? \ldots .	18
	E.3	How to store, share or publish your performed GWAS?	19
	E.4	How to download summary statistics of your results?	20
	E.5	How to download publicly available data?	20

Abstract

Motivation: The rapid growth in genome-wide association studies (GWAS) in plants and animals has brought about the need for a central resource that facilitates i) performing GWAS, ii) accessing data and results of other GWAS, and iii) enabling all users regardless of their background to exploit the latest statistical techniques without having to manage complex software and computing resources.

Results: We present easyGWAS, a web platform that provides methods, tools and dynamic visualizations to perform and analyze GWAS. In addition, easyGWAS makes it simple to reproduce results of others, validate findings, and access larger sample sizes through merging of public datasets. **Availability:** Detailed method and data descriptions as well as tutorials are available in the supplementary materials. easyGWAS is available at http://easygwas.tuebingen.mpg.de/. **Contact:** dominik.grimm@tuebingen.mpg.de

1 Motivation

Genome-wide association studies (GWAS) are an integral tool for discovering the polygenic architecture underlying many complex traits. The recent steady growth of GWAS applications across species (e.g. Arabidopsis thaliana [1, 2, 3], Drosophila melanogaster [4]) has generated a wealth of genotypic and phenotypic data, which makes it possible to search for significant association signals for multiple traits in one species or related traits in different species. Further analysing the shared associations across traits may provide invaluable insights for genetics and evolutionary biology: First, comparing results of GWAS in different species may enable statistical validation of association signals, for instance, to further support findings in human genetics by replicating compatible results in model organisms. Second, one can examine the hypothesis that genetic factors influencing phenotypes can be traced back to related genes and mutations across species. By discovering these common genetic origins of phenotypes, we may gain a deeper understanding of the adaptation of species and possibly of the convergent or parallel evolution of complex traits.

1.1 Difficulties in performing GWAS across traits and species

Obtaining genome-wide association mapping results across multiple studies, traits or even species, however, is still a cumbersome enterprise, which is complicated by three problems: First, several software packages (e.g. PLINK [5]) or species-specific websites (e.g. DGRP [4], Matapax [6], Emma-Server [7]) allow to perform genome-wide association studies on a given dataset. However, they either do not provide genotype and phenotype data at all or only for one single species. Second, existing databases for GWAS results (e.g. GWAS Catalog [8]) focus primarily on human genetics and present summary statistics only for the top scoring loci. This is despite the fact that the most significant genetic loci alone often explain only a small fraction of the heritability of complex traits [9]. Third, databases for human genotypic and phenotypic data such as NCBI dbGaP [10] require a formal application and evaluation before data access can be granted. Data from GWAS

in model organisms and crops is more easily accessible in principle, but can only be obtained from individual websites in a variety of data formats. If one wants to run association studies with identical parameter settings on these datasets, one has to perform tedious data preprocessing and data integration steps first.

1.2 Role of easyGWAS

The field is missing a platform that allows for easy and open access to published genotype and phenotype data from model organisms and crops and is able to perform GWAS on different traits and species. Here, we announce the release of easyGWAS, an interactive and easy-to-use online platform, whose purpose is to fill this gap. easyGWAS is available at http://easygwas.tuebingen.mpg.de/ and enables users to perform GWAS online on their own private or publically available data and to store and publish phenotypic data, meta information on the samples and GWAS results.

1.3 Functionality

Genotypic and phenotypic data from a continuously growing set of published GWA studies are prestored in easyGWAS. Users can either work with these public phenotypes or upload their own phenotypes. Phenotype data uploaded by the user can either be kept private for primary analyses, shared with a restricted set of collaborators, or made publicly available to the community such that other researchers can reuse them in their GWAS analyses. easyGWAS allows to perform univariate association tests in an interactive manner, without the need to manage any computing resources or software. Depending on whether the phenotype is binary or continuous, easyGWAS offers suitable types of mapping algorithms to the user. The results of a completed GWAS can then be visualized in Manhattan plots with gene annotations for the top scoring signals. The example in section 3 and the Supplementary Material includes detailed instructions on how to perform GWAS in easyGWAS.

1.4 Inter-compatibility with statistical genetics software packages

If the user wants to perform an analysis, which is currently only available in existing statistical genetics software packages but not in easyGWAS, the user can export easyGWAS data and store them locally in a file format readable for these software packages. Data export to PLINK, comma-separated files (CSV) and hierarchical data format (HDF5) file format is already available in easyGWAS.

2 Example

In this example of usage, we demonstrate how to perform a GWA study in the plant model organism Arabidopsis thaliana. We use an already published phenotype FLC, which is related to the flowering time of the plant. For the FLC phenotype, RNA was extracted from leaves after four weeks of growth and gene expression levels were determined by northern hybridization quantified relative to

beta–tubulin expression.

Creating a new GWA experiment is divided into several intuitive steps. To be able to follow our instructions, the user should first navigate to the easyGWAS experiment wizard by clicking on **GWAS Center** and then on **Create new GWAS**.

- 1. First, the user selects a species and a dataset. In our example, we select the species Arabidopsis thaliana and the dataset AtPolyDB (call method 75, Horton et al.) and then click on **Continue**.
- 2. In the second step, the user has to select a phenotype. Here the user has the choice to select a published, private/shared or public phenotype. Published phenotypes are accompanied by a peer-reviewed research article, public phenotypes were made public by another user, but need not originate from a publication. Additionally, the user can upload own data (see Supplementary Materials or the online FAQ for detailed tutorials). Here, we select an already published phenotype *FLC* [3]. For this purpose, we select the tab 2.1 Select an existing published phenotype and type into the input field the name of the phenotype *FLC*. Autocompletion will support the user to select the correct phenotype. We proceed by clicking on Continue.
- 3. In this step, users can add additional factors such as principle components or covariates (e.g. environmental factors, gender specific characteristics). In our example, we do not add any additional factors. We click **Continue** in the tab **3.1 No additional factors**.
- 4. Now the user has to select genotypic data. Here, all provided SNPs, specific chromosomes or a region of SNPs can be selected. We select chromosome 1 and 5 in the tab **4.2 Select one or several chromosomes for** *Arabidopsis thaliana* by checking the boxes and click on Select chromosomes.
- 5. To perform a GWAS, we have to select the association method we intend to use. In the algorithms view, one has options to also apply different transformations or filter to the data, such as normalizing the phenotypes. The selection of methods and transformations is dependent on the chosen data. Our web application is analyzing the data on the fly and is enabling only those options that are applicable for the chosen data. Here we keep the default settings using a **Linear Regression** without any transformations. Then we click on **Continue**.
- 6. In the last step, users can check all inputs and can make adjustments if necessary. If everything is correct, the experiment can be submitted to the computation servers. For this purpose, we simply click on **Submit Experiment**.

Finally, the experiment is submitted and all computations are performed in the background. The current view refreshes every 3 seconds. In the meanwhile, users can submit new experiments or browse the data. Nevertheless, this example is finished in around 60 seconds and you will get automatically redirect to the result view to analyze the results. In the result view we provide

dynamic Manhattan plots. Every single SNP can be explored in more detail by moving the mouse over a single point in the plot. On the left we provide a list with the top 10 SNPs and in which gene they are located. In our example the user can see at a glance, that for example the top three SNPs are located in chromosome 5. Additionally, we provide a more detailed SNP annotation view, quantile-quantile plots (QQ-plots) and a phenotype specific view with details about the phenotype (see Supplementary Materials for a detailed description and screenshots). Summary statistics can be downloaded in various formats for further analysis with third party tools.

Additional detailed tutorials (supported by screenshots) about uploading, sharing and downloading data are included in the Supplementary Materials.

3 Conclusion and future plans

easyGWAS is designed to be a dynamically evolving platform with a growing number of functions and prestored datasets. As of now, easyGWAS offers published genotypic and phenotypic data for *Arabidopsis thaliana* [1, 2, 3] and *Drosophila melanogaster* [4] and users can upload their own phenotypic data. easyGWAS enables single-locus mapping with population structure correction for a single trait at a time.

In future versions of easyGWAS, we plan to extend the list of species and to allow users to upload their own genotypic data, while retaining data quality and reliability. Further methods for multi-locus and multi-trait mapping and for automatically retrieving shared association signals across traits will be included.

In summary, we believe that easyGWAS will foster new types of genetic analyses, by providing a convenient framework, which includes data and algorithms for obtaining GWAS results across traits, studies and species.

A Data: Genotypic & phenotypic data and meta information

A.1 Available published data

To easily perform genome-wide association studies (GWAS) across different species a variety of published genotypes and phenotypes are pre-stored in the *easyGWAS* database. As of November 2012, data for *Arabidopsis thaliana* and *Drosophila melanogaster* are available in our public database. For *Arabidopsis thaliana* we included different data sources. The first dataset ['AtPolyDB (call method 75, Horton *et al.*)'] includes 1,307 samples presented by Horton *et al.* in 2012 [1]. Furthermore, we included 107 phenotypes, described and analyzed by Atwell *et al.* [3]. These 107 phenotypes are measured for a subset of these 1,307 samples. The second dataset ['80 genomes data (Cao *et al.*)'] includes 80 samples from the first phase of the 1001 genomes project in *Arabidopsis thaliana* [2]. The genome matrix from the 1001 genomes website¹ was used to retrieve all single nucleotide polymorphisms (SNPs). For this purpose, we excluded all positions with incomplete information and kept all positions with at least one consecutive nucleotide. All SNPs in these *Arabidopsis thaliana* datasets are homozygous ones. Each allele in the SNPs is encoded as described in Table 1.

	major allele	minor allele
major allele	0	1
minor allele	1	2

Table 1: SNP encoding

For the species *Drosophila melanogaster* we integrated a dataset ['Drosophila Genetic Reference Panel (DGRP, Mackay *et al.*)'] with 172 samples [4], sequenced and analyzed by Mackay *et al.*, as well as three phenotypes [4, 11, 12, 13] (six phenotypes, after splitting those into male and female). Missing SNPs in the *Drosophila melanogaster* genome were imputed using the majority allele (different modes of imputation are currently being included into *easyGWAS* and will be available soon).

Additionally we integrated gene annotations for all organisms. This information is used to identify if a SNP is located within a gene or not.

Publicly available genotypes and phenotypes are accompanied by additional meta information such as growth conditions in *Arabidopsis thaliana* or *wolbachia status* in *Drosophila melanogaster*. All datasets were downloaded from their official websites (Table 2).

¹www.1001genomes.org

	Genotypes	https://cynin.gmi.oeaw.ac.at/home/resources/atpolydb
AtPolyDB (call method 75, Horton <i>et al.</i>)	Phenotypes	https://cynin.gmi.oeaw.ac.at/home/resources/atpolydb
	1 nenotypes	http://arabidopsis.gmi.oeaw.ac.at:5000/DisplayResults/
80 genomes data (Cao <i>et al.</i>)	Genotypes	http://1001 genomes.org/data/MPI/MPICao2010/releases/
Arabidopsis thaliana annotations	Annotations	http://www.arabidopsis.org/
Drosophila Genetic Reference	Genotypes	http://dgrp.gnets.ncsu.edu/freeze1/Illumina
Panel (DGRP, Mackay <i>et al.</i>)	Genotypes	$_+_454$ _SNP_genotypes_filtered_for_GWAS/
	Phenotypes	http://dgrp.gnets.ncsu.edu/freeze1/Phenotypes/
Drosophila melanogaster annotations	Annotations	ftp://ftp.flybase.net/releases/FB2008_10/dmel_r5.13/gff/

Table 2: Data sources for all integrated organisms

A.2 How to upload new data

A.2.1 Phenotypic data and meta information

Registered users can upload private data such as phenotypes, covariates or meta information using the easyGWAS wizard (see *Tutorial 5.2*). This data can be used to perform new GWAS or can be shared with collaborators and colleagues. Furthermore, data can made be public to the scientific community. We distinguish between **published** and **public** data. Published data is integrated by the easyGWAS team using data from peer-reviewed publications. However, **public** data was made public by any easyGWAS user. We also provide a contact form for authors who would like to have their published phenotype data integrated into easyGWAS through the easyGWAS team, rather than uploading their data themselves.

A.2.2 New genotypic data

Until now, we provide different datasets for two species. We plan to include datasets from different species to provide a richer selection of genotypic and phenotypic data. To retain quality, we provide an application form which users can use to send a formal data submission application to us. We then will evaluate the request. After successful evaluation, we will provide an upload link and after successful upload and quality inspection of the data, our team will include the data into our database (Figure 1).

A future extension will be a private upload option for small genotypic data sets.

B Integrated methods to perform genome-wide association studies

In general, performing a genome-wide association study is not trivial. There are three main aspects that have to be considered. The first category is data preprocessing. The scientist has to know how to encode, normalize and filter the genotypes, phenotypes and covariates. Second, the scientist has to know which method can be used for which kind of data. There are binary and continuous phenotypes as well as homozygous and heterozygous genotypes. Only specific methods can be applied to a specific type of data. Third, it is crucial to decide whether one should correct for

Steps to submit new Genotype data

New Genotype Data

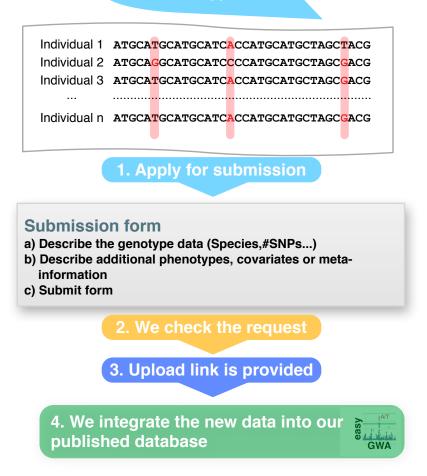


Figure 1: Application process to submit new genotypic data

population stratification or latent confounding factors. Further, some of these methods are hard to parameterize or complicated to set up. One of the strengths of easyGWAS is that it provides several implemented methods and data transformations out of the box. This helps the user to easily perform a genome-wide association study.

B.1 Methods to perform a GWAS

The initial version provides several univariate algorithms, such as linear regression, linear mixed models (EMMAX [14], FaSTLMM [15]) and the Wilcoxon rank-sum test. Linear regression can be used to find single associations between a single SNP and a phenotype. Linear mixed models are used to correct for population structure, family structure and cryptic relatedness at the same

time. To all these methods one can add covariates, such as principle components (PCs), environmental factors or gender specific characteristics. Additionally, the Wilcoxon rank-sum test can be used for homozygous genotypes. These methods are state-of-the-art, more methods will be added continuously. New methods for multi-marker discovery such as two-locus search using graphical computing units (GPUs) and multi-trait discovery will be added in the near future.

B.2 Transformations to standardize data

To transform phenotypic and genotypic data we added several methods. Genotypes can be standardized, one can zero-mean the data and/or divide by unit variance. Phenotypes can be transformed in the same way. Additionally, phenotypes can be log-transformed, square root and box-cox transformed. Figure 2 illustrates a scheme of all options.

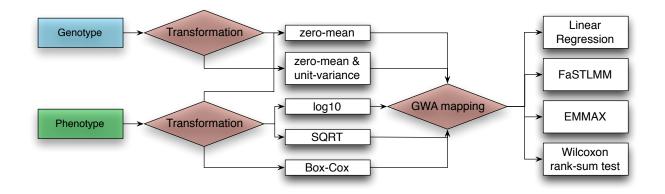


Figure 2: Scheme off all possible transformations and GWA mapping methods

C The web application interface

The web application contains three main parts. There is one view to plan, perform and store GWAS, a second view to browse and analyze the data and a third view to download available datasets. In the following, we will describe all sections in detail.

C.1 The easyGWAS wizard and experiment history

The first section contains all necessary tools to plan, perform and analyze whole genome-wide association studies. Here registered users can use a step-by-step procedure (software wizard) to easily create new experiments (see Tutorial 5.1, Figure 3a). The wizard is divided into several steps. First the user has to choose an available species and dataset. In the next step a single phenotype can be selected. Here it is possible to select an already published, private or public phenotype. Additionally, one can upload an own phenotype. We distinguish between **published** and **public** phenotypes. **Published** phenotypes are already published, whereas **public** phenotypes

a)	CREATE EXPERIMENTS	GWA Wizard 1. Species / 2. Phe	notypes / 3. Additional Factors / 4. Genotypes / 5. Algorithm / 6. Summary	
b)	Crutorial Experiment overview	1. Select a specie	es and a dataset for your GWA experiment	
c)	 O My temporary history ▲ My experiments C Shared experiments C Public experiments 	Select a species: Select a dataset:	Arabidopsis thaliana APOlyDB (call method 75, Horton et al.) \$	
			Continue	

Figure 3: Screenshot from the GWA experiment view. Creating new GWA experiments and sharing the results with collaborators or the scientific community

are phenotypes uploaded by a specific user and made public to the community. After chosen a phenotype it is possible to add additional factors to the experiment, such as principle components or one or several covariates. Covariates are meta information such as environmental factors or gender specific characteristics. To proceed the user has to select genotypic data. To do so, it is possible to select all genotypic data, meaning all available SNPs. Furthermore, specific chromosomes or a range of SNPs can be selected. The last step provides different algorithms, standardizations and filters. The selection of methods is based on the selected genotypic and phenotypic data in the previous steps. The summary view in the end provides all user specific selections and offers the user to adjust settings or to submit the experiment to the computation workers.

Each experiment performed is saved in a temporary experiment history (Figure 3b). Here all experiments are stored for primary analysis for at least 48h. To keep interesting findings, users can store experiments permanently in their private profile. To simplify scientific exchange, all experiments can be shared via the web application with collaborators and colleagues (Figure 3c). Sharing experiments and data can prevent laborious extracting of data and findings. Furthermore, data and experiments can be made public to the scientific community. All summary statistics can be downloaded to further analyze the data using third party tools.

Performing GWAS can be time consuming. Due to an advanced technology (see *Web application infrastructure*) the user can continue working using the web application while the experiment is computed in the background at the same time. Automated email notifications are send out as soon the computations are done. Additionally, the user can track the status of all experiments through the temporary history (Figure 4).

To examine individual experiments each experiment has an interactive results page. Figure 5 shows a screenshot of the result page. The view is divided into two parts. The left part provides general information. Here a short summary table informs about all settings made by the user, e.g. which species, dataset and parameters were selected (Figure 5a). At a glance the user can see the top 10 SNP annotations with the smallest p-values (Figure 5b). Dynamic Manhattan-plots for all chromosomes are rendered in the right half (Figure 5c). Each SNP within the Manhattan plot is interactive, meaning that the user is able to inspect single SNPs getting live information like the corresponding p-value or in which gene the SNP is located. The green line in each Manhattan-

REATE EXPERIMENTS									
Overview									
Create new GWA	Temporary e	xperiments are available for 4	8h! To store them	permanently ple	ase save	the experimen	t using the sa	ve button.	
Tutorial									
XPERIMENT OVERVIEW	Experiment	Experiment id	Date	Species	SNPs	Phenotypes	Algorithm	Running	Operation
My temporary history	name								
My experiments	Experiment 0	299ee52c-182e-49b7-b6a5-	Nov. 11, 2012,	Arabidopsis	all	Leaf serr 16	Linear	0	
Shared experiments		f10699ce7b3e	05:08 PM	thaliana			Regression		

Figure 4: Screenshot from temporary experiment history. The red highlighted row indicates that the experiment is still computing.



Figure 5: Screenshot showing the result page of an experiment.

plot is the Bonferroni threshold. The alpha significance level can be adjusted using the plotting options (Figure 5d). The strength of population structure confounding can be easily explored with

Q-Q plots (Figure 5e) and the genomic control λ . To see the actual distribution of a phenotype, the *Phenotype Explorer* shows histograms for transformed and non-transformed phenotypes and computes a Shapiro-Wilk test to test the null hypothesis that the data was drawn from a normal distribution (Figure 5f). To examine if SNPs of interest are located within genes, significant loci are summarized in a gene-annotation view (Figure 5g).

C.2 The data center

The second main section of the web application is the **Data Center**. Here, the user can browse available data, such as samples, phenotypes and covariates. Detailed information can be accessed for each data entry, such as meta and/or geographical information (Figure 6). Associated publications are provided for all published entries. The **Data Center** contains two main views. One for **published** and **public** data and a second for all user specific **private/shared** data. Private data is only visible to the owner of the data. Note that privately shared data belongs to the owner, meaning that only the owner has the permission to delete or modify shared data.

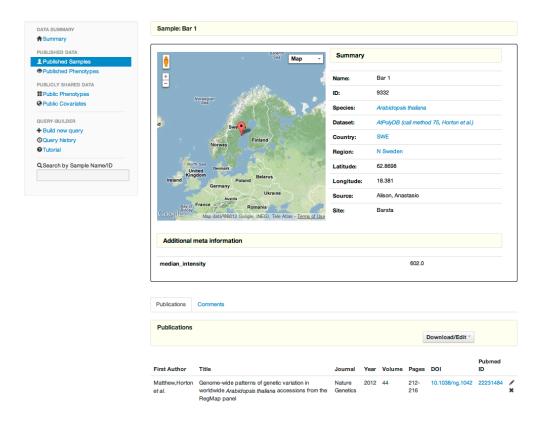


Figure 6: Data center with detailed information about sample, phenotypes and covariates

C.3 The download center

The third section provides additional download options. Here, whole datasets (genotypic and phenotypic data) for all integrated species can be downloaded in different file formats. At the moment we provide the following formats: PLINK[5], comma-separated files (CSV) and hierarchical data format 5 (HDF5)².

D The web application backend

The backend of the web application is completely written in Django³ a web framework for Python. For the web design, we used the Cascading Style Sheets (CSS), provided by Twitter Bootstrap⁴. To

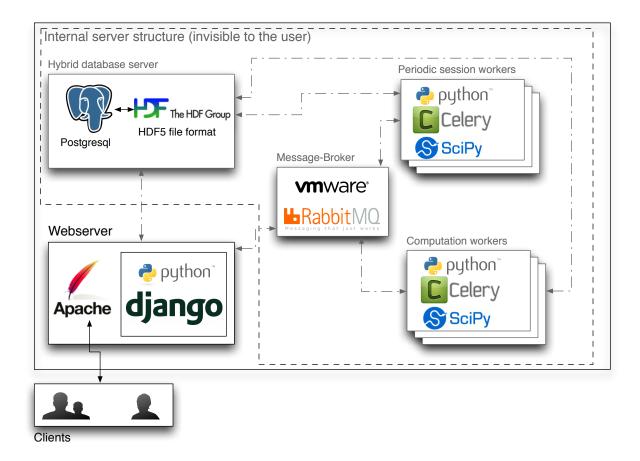


Figure 7: Scheme of the web application infrastructure

handle the huge amount of SNP data we developed a hybrid database model using a PostgreSQL database and the $HDF5^2$ file format. Here all SNPs, associated positions and chromosome indices as well as all phenotypes and covariates are stored in the HDF5 file. Additional phenotype, sample,

²http://www.hdfgroup.org/HDF5/

³https://www.djangoproject.com

 $^{^{4}}$ http://twitter.github.com/bootstrap/

covariate and meta information are stored in the PostgreSQL database. HDF5 files are highly optimized to handle huge files and can be accessed fast and easily. As GWA mappings are resource consuming all computations are distributed to different computation servers (workers). To schedule different tasks smartly we are using a message broker (RabbitMQ⁵). This broker distributes the different tasks to single workers (Figure 7). The backend is well designed to easily extend the functionality of *easyGWAS* using additional novel methods. New species and datasets can be integrated within hours, depending on the size of the data. If more computational power is needed, new workers can be added dynamically.

E Tutorials

In this section we provide various tutorials on how to use easyGWAS. We demonstrate how to actually perform a genome-wide association mapping, how to upload own private phenotypes and how to share or make them public for collaborators and/or the scientific community. Furthermore, we show how to download summary statistics of GWA experiments and published genotype and phenotype data. Screenshots are attached to all important steps.

E.1 How to perform a GWAS easily?

In this tutorial we demonstrate how to easily perform a GWA study using Arabidopsis thaliana and already published phenotype FLC (flowering time related phenotype).

- 1. If not already done: Create a new $easyGW\!AS$ account and log in.
- 2. Navigate to the easyGWAS wizard Menu: GWA-Experiments \rightarrow Create new GWA
- 3. First, select a species and a dataset. Here we choose the species Arabidopsis thaliana and the dataset AtPolyDB (call method 75, Horton et al.)[1] and click **Continue**.

CREATE EXPERIMENTS	GWA Wizard
+ Create new GWA Tutorial	1. Species / 2. Phenotypes / 3. Additional Factors / 4. Genotypes / 5. Algorithm / 6. Summary
EXPERIMENT OVERVIEW O My temporary history My experiments C Shared experiments C Public experiments	1. Select a species and a dataset for your GWA experiment Select a species: Arabidopsis thaliana Select a dataset: AlPolyDB (call method 75, Horton et al.) Continue

4. Select a phenotype. Here the user has the choice to select published, private/shared and public phenotypes. Additionally the user can upload his own data (see Tutorial 5.2). Here

⁵http://www.rabbitmq.com

we select a published phenotype FLC[3]. For this purpose, select the tab **2.1 Select an** existing published phenotype and type into the input field the name of the phenotype (FLC[3]). Auto-completion will help you to select the correct one. Click Continue.

Control Contro Control Control Control Control Control Contro	CREATE EXPERIMENTS	GWA Wizard
Crutorial Excensioner overview OMy temporary history My experiments Select appendent overview Select binary phenotype Select continuous phenotype Select continuous phenotype Select by name Select binary phenotype Select continuous phenotype Select appendent overview Select appenden	A Overview	
EXPERIMENT OVERVIEW OMy temporary history My experiments Charact experiments Charact experiments Select by name Select by name Select by name Select by name Select a phenotype by name among all phenotypes (autocompletion support): FL ArboyDB (call method 75, Horton et al.) 2.2 Upload new phenotypes 2.3 Select a private or shared phenotype	+ Create new GWA	1. Species / 2. Phenotypes / 3. Additional Factors / 4. Genotypes / 5. Algorithm / 6. Summary
Contract of the second	Tutorial	
EMy experiments Charded experiments Charded experiments Select by name	EXPERIMENT OVERVIEW	2. Select phenotypes for your GWA experiment
C Shared experiments C Public experiments C Public experiments C Public experiments C Public experiments C Select by name Select binary phenotype Select continuous phenotype Select continuous phenotype Continue C Continue C Continue C Continue C C C C C C C C C C C C C C C C C C C	OMy temporary history	
© Public experiments Select binary phenotype Select continuous phenotype LXPERIMENT OVERVIEW Select a phenotype by name among all phenotypes (autocompletion support): FL FL FLC FL FLC 22 Upload new phenotypes 2.3 Select a private or shared phenotype 2.3 Select a private or shared phenotype	My experiments	2.1 Select an existing published phenotype
EXPERIMENT OVERVIEW 1. Selected species: Arabid/opsis trailana APolyDB (call method 75, Horton et al.) 2.2 Upload new phenotypes 2.3 Select a private or shared phenotype	C Shared experiments	
EXPERIMENT OVERVIEW I. Selected species: Arb/dy/Se (call method 75, Horton ef al.) 2.2 Upload new phenotypes 2.3 Select a private or shared phenotype	Public experiments	Select by name Select binary phenotype Select continuous phenotype
2.3 Select a private or shared phenotype	1. Selected species: Arabidopsis thaliana AtPolyDB (call method 75,	FL Continue
		2.2 Upload new phenotypes
		2.3 Select a private or shared phenotype
2.4 Select a <i>public</i> phenotype		2.4 Select a public phenotype

5. In the following step you can add additional factors such as principle components or covariates (e.g. environmental factors, gender specific characteristics). In this tutorial we do not add any additional factors. Click **Continue** in the tab **3.1 No additional factors**.

CREATE EXPERIMENTS	GWA Wizard
Overview	
+ Create new GWA	1. Species / 2. Phenotypes / 3. Additional Factors / 4. Genotypes / 5. Algorithm / 6. Summary
Tutorial	
EXPERIMENT OVERVIEW	3. Select additional factors such as principle components or covariates
OMy temporary history	
My experiments	3.1 No additional factors
C Shared experiments	
Public experiments	Continue if you do not like to add additional factors.
EXPERIMENT OVERVIEW	3.2 Covariate - Principle component analysis (PCA)
1. Selected species:	3.3 Covariate - Upload new covariates
Arabidopsis thaliana AtPolyDB (call method 75,	
Horton et al.)	
	Back
2. Selected phenotype:	
FLC	

6. Now we have to select genotypic data. You have to choose if you like to use all provided SNPs, specific chromosomes or a region of SNPs. For this tutorial we select chromosome 1 and 5 in the tab **4.2 Select one or several chromosomes for** *Arabidopsis thaliana* by checking the boxes. Click **Select chromosomes**.

CREATE EXPERIMENTS	GWA Wizard
A Overview	
+ Create new GWA	1, Species / 2, Phenotypes / 3, Additional Factors / 4, Genotypes / 5, Algorithm / 6, Summary
Tutorial	
EXPERIMENT OVERVIEW	4. To select SNPs for your GWA experiment choose one of the following options
OMy temporary history	
My experiments	4.1 Select all SNPs for Arabidopsis thaliana
C Shared experiments	4.2 Select one or several chromosomes for Arabidopsis thaliana
Public experiments	4.2 detections of several circlinosomes to Arabicopas tranana
	Select one or several chromosomes from Arabidopsis thaliana
EXPERIMENT OVERVIEW	Chromosome Chr1
	Chromosome Chr2
1. Selected species: Arabidopsis thaliana	Chromosome Chr3
AtPolyDB (call method 75,	
Horton et al.)	Chromosome Chr4
	Chromosome Chr5 Select chromosomes
2. Selected phenotype:	
FLC	4.3 Select a subset of SNPs for Arabidopsis thaliana
3. Additional Factors:	
None	Back

7. To perform a GWAS we have to select a method we intend to use. In the algorithms view one has options to also apply different transformations or filter to the data. The selection of methods and transformations is dependent on the chosen data. Our web application is analyzing the data on the fly and is enabling only those options that are applicable for your data. Here we keep the default settings using a **Linear Regression** without any transformations. Click **Continue**.

Overview		
Create new GWA Tutorial	1. Species / 2. Phenotypes / 3. Additional Factors / 4. Genotypes / 5. Algori	ithm / 6. Summary
XPERIMENT OVERVIEW	5. Select an algorithm and additional parameters for your GWA experi	iment
My temporary history		
My experiments	Provided choice of algorithms and standardizations depends on the se	lected phenotype. More information in the FAQ.
Shared experiments		
Public experiments	5.1 Algorithm	
EXPERIMENT OVERVIEW	Choose an algorithm from the list	Linear Regression \$
1. Selected species: Arabidopsis thaliana AtPolyDB (call method 75, Horton et al.)	Algorithm Description	
Arabidopsis thaliana	Algorithm Description Linear regression is used to test if a arbitrary phenotype is associated with a marker/SNP is tested individually. Additional factors/covariates can be adde specific covariates or principle components (PCs).	
Arabidopsis thalana AtPolyDB (call method 75, Horton <i>et al.</i>) 2. Selected phenotype: FLC	Linear regression is used to test if a arbitrary phenotype is associated with a marker/SNP is tested individually. Additional factors/covariates can be adde	
Arabidopsis thaliana AtPolyDB (call method 75, Horton et al.) 2. Selected phenotype:	Linear regression is used to test if a arbitrary phenotype is associated with a marker/SNP is tested individually. Additional factors/covariates can be adde specific covariates or principle components (PCs).	
Arabidopsis malana AIPohyon B (call method 75, Horton et al.) 2. Selected phenotype: FLC 3. Additional Factors:	Linear regression is used to test if a arbitrary phenotype is associated with a marker/SNP is tested individually. Additional factors/covariates can be adde specific covariates or principle components (PCs).	d to the model such as environmental factors, gender
Arabidopois finalana AtPolyDB (call method 75, Horton et al.) 2. Selected phenotype: FLC 3. Additional Factora: None 4. Selected SNPs:	Linear regression is used to test if a arbitrary phenotype is associated with a marker/SNP is tested individually. Additional factors/covariates can be adde specific covariates or principle components (PCs).	d to the model such as environmental factors, gender
Arabidopois finalana AtPolyDB (call method 75, Horton et al.) 2. Selected phenotype: FLC 3. Additional Factora: None 4. Selected SNPs:	Linear regression is used to test if a arbitrary phenotype is associated with a marker/SNP is tested individually. Additional factors/covariates can be adde specific covariates or principle components (PCs). 5.2 Filter Choose threshold for minimum allele frequency (mAF) 5.3 Standardization	d to the model such as environmental factors, gender
Arabidopois finalana AtPolyDB (call method 75, Horton et al.) 2. Selected phenotype: FLC 3. Additional Factora: None 4. Selected SNPs:	Linear regression is used to test if a arbitrary phenotype is associated with a marker/SNP is tested individually. Additional factors/covariates can be adde specific covariates or principle components (PCs). 5.2 Filter Choose threshold for minimum allele frequency (mAF) 5.3 Standardization Choose a standardization method for your genotypes	d to the model such as environmental factors, gender

8. In the last step you can check all your inputs again and do necessary adjustments. If everything is correct, you can submit your experiment to the computation servers. For this purpose, simply click **Submit Experiment**.

9. Finally your experiment is submitted. All computations are running in the background. The current view gets refreshed every 3 seconds. In the meanwhile you could submit new experiments or browse the data. Nevertheless, this experiment is finished in around 60 seconds and you will get automatically redirect to the result view.

Overview	
+ Create new GWA	Providence and wellet all
Tutorial	Experiment submitted!
EXPERIMENT OVERVIEW	If the computations are ready you will get notified via a push notification or you directly got redirected to the results-page. If you leave the web-application you get notified via e-mail.
OMy temporary history	ieave me web-application you get notified via e-mail.
My experiments	Your gwa-experiment is saved in your temporary history.
Shared experiments	
Public experiments	Unique Experiment-ID: db0770fc-f06c-4d78-967a-3847bb57c8ae
	0
	Experiment is running 15.0 seconds

E.2 How to upload an own phenotype and perform a GWAS on it?

Here we show how to easily upload new phenotypic data and how to perform a GWAS with it.

- 1. Navigate to the easyGWAS wizard Menu: GWA-Experiments \rightarrow Create new GWA
- 2. First, select a species and a dataset. Here we choose the species Arabidopsis thaliana and the dataset AtPolyDB (call method 75, Horton et al.)[1] and click **Continue**.

CREATE EXPERIMENTS	GWA Wizard						
+ Create new GWA	1. Species / 2. Phenotypes / 3. Additional Factors / 4. Genotypes / 5. Algorithm / 6. Summary						
EXPERIMENT OVERVIEW O My temporary history My experiments C Shared experiments C Public experiments	1. Select a species and a dataset for your GWA experiment Select a species: Arabidopsis thaliana Select a dataset: AtPolyDB (call method 75, Horton et al.)						

3. Now navigate to **2.2 Upload new phenotypic data** and download the linked **demo** file. Click on **Choose File** and upload the demo file. Click on **Continue**.

GWA Wizard						
1. Species / 2. Phenotypes / 3. Additional Factors / 4. Genotypes / 5. Algorithm / 6. Summary						
2. Select phenotypes for your GWA experiment						
2.1 Select an existing published phenotype						
2.2 Upload new phenotypic data						
Upload a new phenotype. To check the file-format or get additional informations please see the FAQ.						
To browse all available accessions click here: Browse accessions						
Here you can download an example file. Download demo file (Arabidopsis thalians)						
Upload file with new phenotype:						
Upload file with new phenotype:						
Upload file with new phenotype: Choose File no file selected						
Upload file with new phenotype: Choose File no file selected						
Upload file with new phenotype: Choose File no file selected Continue						
Upload file with new phenotype: Choose File no file selected Continue 2.3 Select a private or shared phenotype						
Upload file with new phenotype: Choose File no file selected Continue 2.3 Select a private or shared phenotype						

4. Finally proceed like in Tutorial 5.1. You may try different algorithms and transformations.

E.3 How to store, share or publish your performed GWAS?

To keep interesting findings and experiments users can save their results permanently using the saving functionality in **My temporary history**. Here you can rename your experiment and save it in your experiment history **My experiments**.


Temporary 0	WA-Experiment history	Saved GWA-Exp	periments S	hared GWA-Expe	eriments			
Temporary experiments are available for 48h! To store them permanently please save the experiment using the save button.								
Experiment name	Experiment id	Date	Species	SNPs	Phenotypes	Algorithm	Running	Operations
Experiment 0	738ae045-da54-4a05- a6f0-409563e7d931	Nov. 11, 2012, 06:32 PM	Arabidopsis thaliana	chromosome	FLC	Linear Regression	Done	×∕⊚ ಔ≞
Experiment 1	fbf6a052-efaf-4725-a79d- 0d9b4cf66c91	Nov. 11, 2012, 06:34 PM	Arabidopsis thaliana	chromosome	FLC	Linear Regression	Done	ave experiment

To simplify data exchange between colleagues and the scientific community one has the possibility to share and publish saved experiments. For this purpose, open the category **My experiments** and click on either **Share Experiment** or **Publish experiment**. Note that if you like to publish an experiment which was performed on private phenotypes and/or covariates, one has to publish all dependent private data. Please provide meaningful and useful names.

Make your experiment public to the scientific community							
WATCH OUT! Published experiments cannot be deleted anymore. Please provide a meaningful name!							
Private phenotypes and private covariates are linked to this experiment. If you like to publish this experiment you also have to publish all associated phenotypes and covariates. Please provide meaningful names!							
Experiment name:	Experiment 1						
Experiment id:	f2888e0b-f1f7-4da8-9c24-118376536d9c	- Drouido mooninaful nomoo					
Covariate (id: AT_C_0):	TestCovariate	Provide meaningful names					
Covariate (id: AT_C_1):	TestCovariate2						
Phenotype (id: AT_P_284):	TestPhenotype LD	renotype LD					
		Publish					

E.4 How to download summary statistics of your results?

For each experiment summary statistics can be downloaded using the download assistant at the GWAS result page. Click on **Download Summary Statistics** and choose one of your preferred formats. Right now there are choices for comma-separated files (CSV) and hierarchical data format 5 (HDF5²). The summary statistic files contain p-values for each loci and chromosome. The HDF5 file has additional information on which samples were used.

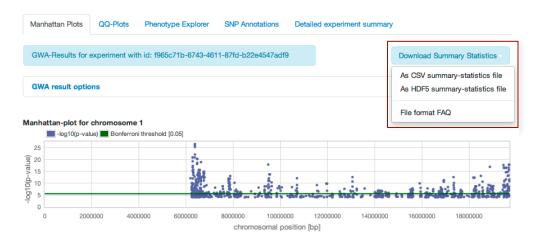

E.5 How to download publicly available data?

To download whole genotype and phenotype data you can use the **Download Center** in the main menu. Here you can download all available data for each species and dataset in various formats, supporting PLINK[5], CSV and HDF5 format.

References

- Horton, M. W., Hancock, A. M., Huang, Y. S., Toomajian, C., Atwell, S., Auton, A., Muliyati, N. W., Platt, A., Sperone, F. G., Vilhjalmsson, B. J., Nordborg, M., Borevitz, J. O., and Bergelson, J. Nat Genet 44(2), 212–216 February (2012).
- [2] Cao, J., Schneeberger, K., Ossowski, S., Gnther, T., Bender, S., Fitz, J., Koenig, D., Lanz, C., Stegle, O., Lippert, C., Wang, X., Ott, F., Mller, J., Alonso-Blanco, C., Borgwardt, K., Schmid, K. J., and Weigel, D. *Nature Genetics* 43(10), 956–963 October (2011). PMID: 21874002.
- [3] Atwell, S., Huang, Y. S., Vilhjlmsson, B. J., Willems, G., Horton, M., Li, Y., Meng, D., Platt, A., Tarone, A. M., Hu, T. T., Jiang, R., Muliyati, N. W., Zhang, X., Amer, M. A., Baxter, I., Brachi, B., Chory, J., Dean, C., Debieu, M., Meaux, J. d., Ecker, J. R., Faure, N., Kniskern, J. M., Jones, J. D. G., Michael, T., Nemri, A., Roux, F., Salt, D. E., Tang, C., Todesco, M., Traw, M. B., Weigel, D., Marjoram, P., Borevitz, J. O., Bergelson, J., and Nordborg, M. *Nature* 465(7298), 627–631 March (2010).
- [4] Mackay, T. F. C., Richards, S., Stone, E. A., Barbadilla, A., Ayroles, J. F., Zhu, D., Casillas, S., Han, Y., Magwire, M. M., Cridland, J. M., Richardson, M. F., Anholt, R. R. H., Barrn, M., Bess, C., Blankenburg, K. P., Carbone, M. A., Castellano, D., Chaboub, L., Duncan, L., Harris, Z., Javaid, M., Jayaseelan, J. C., Jhangiani, S. N., Jordan, K. W., Lara, F., Lawrence, F., Lee, S. L., Librado, P., Linheiro, R. S., Lyman, R. F., Mackey, A. J., Munidasa, M., Muzny, D. M., Nazareth, L., Newsham, I., Perales, L., Pu, L., Qu, C., Rmia, M., Reid, J. G., Rollmann, S. M., Rozas, J., Saada, N., Turlapati, L., Worley, K. C., Wu, Y., Yamamoto, A., Zhu, Y., Bergman, C. M., Thornton, K. R., Mittelman, D., and Gibbs, R. A. Nature 482(7384), 173–178 February (2012). PMID: 22318601.
- [5] Purcell, S., Neale, B., Todd-Brown, K., Thomas, L., Ferreira, M. A. R., Bender, D., Maller, J., Sklar, P., de Bakker, P. I. W., Daly, M. J., and Sham, P. C. American journal of human genetics 81(3), 559–575 September (2007). PMID: 17701901.
- [6] Childs, L. H., Lisec, J., and Walther, D. Plant Physiology 158(4), 1534–1541 April (2012).
- [7] Bennett, B. J., Farber, C. R., Orozco, L., Min Kang, H., Ghazalpour, A., Siemers, N., Neubauer, M., Neuhaus, I., Yordanova, R., Guan, B., Truong, A., Yang, W.-p., He, A., Kayne, P., Gargalovic, P., Kirchgessner, T., Pan, C., Castellani, L. W., Kostem, E., Furlotte, N., Drake, T. A., Eskin, E., and Lusis, A. J. *Genome Research* 20(2), 281–290 January (2010).
- [8] Hindorff, L. A., Sethupathy, P., Junkins, H. A., Ramos, E. M., Mehta, J. P., Collins, F. S., and Manolio, T. A. Proceedings of the National Academy of Sciences of the United States of America 106(23), 9362–9367 June (2009). PMID: 19474294.

- [9] Manolio, T. A., Collins, F. S., Cox, N. J., Goldstein, D. B., Hindorff, L. A., Hunter, D. J., McCarthy, M. I., Ramos, E. M., Cardon, L. R., Chakravarti, A., Cho, J. H., Guttmacher, A. E., Kong, A., Kruglyak, L., Mardis, E., Rotimi, C. N., Slatkin, M., Valle, D., Whittemore, A. S., Boehnke, M., Clark, A. G., Eichler, E. E., Gibson, G., Haines, J. L., Mackay, T. F. C., McCarroll, S. A., and Visscher, P. M. *Nature* 461(7265), 747–753 October (2009).
- [10] Mailman, M. D., Feolo, M., Jin, Y., Kimura, M., Tryka, K., Bagoutdinov, R., Hao, L., Kiang, A., Paschall, J., Phan, L., Popova, N., Pretel, S., Ziyabari, L., Shao, Y., Wang, Z. Y., Sirotkin, K., Ward, M., Kholodov, M., Zbicz, K., Beck, J., Kimelman, M., Shevelev, S., Preuss, D., Yaschenko, E., Graeff, A., Ostell, J., and Sherry, S. T. *Nature genetics* **39**(10), 1181–1186 October (2007). PMID: 17898773 PMCID: PMC2031016.
- [11] Harbison, S. T., Yamamoto, A. H., Fanara, J. J., Norga, K. K., and Mackay, T. F. C. Genetics 166(4), 1807 –1823 April (2004).
- [12] Morgan, T. J. and Mackay, T. F. C. *Heredity* **96**(3), 232–242 March (2006). PMID: 16404413.
- [13] Jordan, K. W., Carbone, M. A., Yamamoto, A., Morgan, T. J., and Mackay, T. F. Genome Biology 8(8), R172 August (2007).
- [14] Kang, H. M., Sul, J. H., Service, S. K., Zaitlen, N. A., Kong, S.-y., Freimer, N. B., Sabatti, C., and Eskin, E. Nat Genet 42(4), 348–354 April (2010).
- [15] Lippert, C., Listgarten, J., Liu, Y., Kadie, C. M., Davidson, R. I., and Heckerman, D. Nat Meth 8(10), 833–835 October (2011).