

Enriching and Characterizing T-Cell Repertoires from 3' Barcoded Single-Cell Whole Transcriptome Amplification Products

Authors: Tasneem Jivanjee^{1,2,3,4,#}, Samira Ibrahim^{1,2,3,4,#}, Sarah K. Nyquist^{1,2,3,4,5,6}, G. James Gatter^{1,2,3,4}, Joshua D. Bromley^{1,2,3,4,7}, Swati Jaiswal⁸, Bonnie Berger^{6,9}, Samuel M. Behar⁸, J. Christopher Love^{2,3,4,10,*}, Alex K. Shalek^{1,2,3,4,5*}

[#]These first authors contributed equally to this work.

^{*}These senior authors contributed equally to this work.

Affiliations:

1. Institute for Medical Engineering & Science and Department of Chemistry, Massachusetts Institute of Technology, Cambridge, MA
2. Koch Institute for Integrative Cancer Research, Massachusetts Institute of Technology, Cambridge, MA
3. Ragon Institute of MGH, MIT, and Harvard, Cambridge, MA
4. Broad Institute of MIT and Harvard, Cambridge, MA, Massachusetts Institute of Technology, Cambridge, MA
5. Program in Computational and Systems Biology, Massachusetts Institute of Technology, Cambridge, MA, USA
6. Computer Science and Artificial Intelligence Laboratory, Massachusetts Institute of Technology, Cambridge, MA, USA.
7. Microbiology Graduate Program, Massachusetts Institute of Technology, Cambridge, MA 02139, USA
8. Department of Microbiology and Physiological Systems, University of Massachusetts Medical School, Worcester, Massachusetts, USA

9. Department of Mathematics, Massachusetts Institute of Technology, Cambridge, MA

10. Department of Chemical Engineering, Massachusetts Institute of Technology,
Cambridge, MA

Corresponding authors: **J. Christopher Love:** clove@mit.edu; **Alex K. Shalek:**

shalek@mit.edu

Abstract

Antigen-specific T cells play an essential role in immunoregulation and many diseases such as cancer. Characterizing the T cell receptor (TCR) sequences that encode T cell specificity is critical for elucidating the antigenic determinants of immunological diseases and designing therapeutic remedies. However, methods of obtaining single-cell TCR sequencing data are labor and cost intensive, typically requiring both cell sorting and full length single-cell RNA-sequencing (scRNA-seq). New high-throughput 3' cell-barcoding scRNA-seq methods can simplify and scale this process; however, they do not routinely capture TCR sequences during library preparation and sequencing. While 5' cell-barcoding scRNA-seq methods can be used to examine TCR repertoire at single-cell resolution, doing so requires specialized reagents which cannot be applied to samples previously processed using 3' cell-barcoding methods.

Here, we outline a method for sequencing *TCR α* and *TCR β* transcripts from samples already processed using 3' cell-barcoding scRNA-seq platforms, ensuring TCR recovery at a single-cell resolution. In short, a fraction of the 3' barcoded whole transcriptome amplification (WTA) product typically used to generate a massively parallel 3' scRNA-seq library is enriched for TCR transcripts using biotinylated probes, and further amplified using the same universal primer sequence from WTA. Primer extension using TCR V-region primers and targeted PCR amplification using a second universal primer results in a 3' barcoded single-cell *CDR3*-enriched library that can be sequenced with custom sequencing primers. Coupled with 3' scRNA-seq of the same WTA, this method enables simultaneous analysis of single-cell transcriptomes and TCR sequences which can help interpret inherent heterogeneity among antigen-specific T cells

and salient disease biology. The method presented here can also be adapted readily to enrich and sequence other transcripts of interest from both 3' and 5' barcoded scRNA-seq WTA libraries.

Key words: Single-cell RNA-sequencing, scRNA-seq, Seq-Well S³, T-cell receptor repertoire profiling, Gene expression, Targeted Enrichment, single-cell TCR sequencing

1. Introduction

Antigen-specific T cells play a critical role in immunoregulation and various diseases such as cancer and autoimmune disorders (1–4). Each T cell possesses a genetically encoded T cell receptor (TCR); at any given moment, the repertoire of all TCRs and the T cell states to which each is tied define which specific antigens a host can identify and how it will respond to them, respectively, informing key aspects of adaptive immunity (5). Critically, a host's overall TCR repertoire is not fixed, but rather can change in response to internal and external factors, such as spontaneous mutations or infection, respectively (6). Thus, sequencing this region is essential for characterizing TCR diversity. Yet, the diversity of these TCRs can also present a challenge when it comes to precisely analyzing the repertoire and its relation to cellular phenotypes (7). To improve the accuracy of identifying key T cell subpopulations, methods that can link genetically encoded receptor identity and phenotype to adaptive immune responses are needed (8).

The TCR itself is a heterodimer of two chains ($\alpha\beta$ or $\gamma\delta$). The uniqueness of each T cell's TCR is the result of recombination events occurring between the genes encoding for the TCR α and β (or γ and δ) variable (V), diversity (D), and joining (J) regions. For the TCR α chain, recombination results in one of myriad VJ sequences; for the β , it yields a VDJ. The unique part of each α and β

TCR chain is contained within its complementarity-determining region 3 (CDR3). This region includes the V(D)J recombination junction as well as the part of the receptor that binds to major histocompatibility complex (MHC) molecules on antigen-presenting cells (5).

To investigate the sequences of the TCR α and β chains belonging to specific T cells, a few methods have emerged. One approach is bulk genome sequencing of the entire TCR repertoire. Although this strategy (and related ones) enables comprehensive population-level TCR analysis, it mainly provides insight into overall TCR α and TCR β usage rather than revealing distinct single-cell TCR α -TCR β pairs (though this can be computationally inferred); similarly, it does not allow exploration of gene expression at the single-cell level, hindering understanding of the relationship between TCR clonotype (a cell's specific α , β pair) and biological function (6).

Another analysis method uses tetramers to tag antigen-specific T cells. Here, fluorescently- (or oligonucleotide-) labeled MHC-peptide tetramers are used to mark antigen-specific T cells for subsequent analysis (9). This approach can be combined with cell sorting methods to identify the clonotypes that bind to the specific tetramer. Specifically, tetramer-labeled cells are single-cell sorted and then each cell's $\alpha\beta$ TCR pair is amplified by targeted PCR and determined by sequencing (9, 10). To link clonotype to phenotype, this initial targeted PCR can be replaced by whole transcriptome amplification (WTA); after, the sequences of the *TCR α* and *TCR β* can be reconstructed computationally from full length single-cell RNA-sequencing (scRNA-seq) data or probed directly by targeted library generation and sequencing (11). However, these plate-based methods have limited throughput due to labor, time and cost constraints (6, 12, 13).

To scale phenotypic profiling, newer high-throughput scRNA-seq methods have been developed that rely on early cell barcoding. Early tagging of each cell's mRNAs with unique nucleic acid sequences allows many cells to be processed at once while still enabling computational grouping of the transcripts belonging to each cell (14). While these methods may be useful for high-throughput single-cell gene expression profiling, most 3' barcoded scRNA-seq methods do not provide information on the variable region of the α and β TCR genes—particularly the complementarity-determining region 3 (*CDR3*) regions which lie closer to the 5' end of the TCR transcripts and are effectively lost in most 3' library preparation methods which generate short fragments (~300-500 bp) for sequencing (Fig. 1). In other words, 3' scRNA-seq methods (such as Drop-seq, Seq-Well S³, 10x 3', and inDrops, or similar spatial methods such as Slide-Seq) typically do not return fragments for sequencing with both *CDR3* sequences and the 3' cell barcode and unique molecular identifier (UMI). This information is necessary to determine the identify of a cell's α and β TCR sequences and to link them to a specific cell identifier, respectively, so that they may be tied to a specific phenotypic state (as determined by the other transcripts that the cell expresses) (15–17). While targeted amplification using a primer against the TCR constant region works for 5' cell barcoded WTA products, such as those created using the 10x 5' v2 kit, it requires special reagents (a TCR sequencing kit). Furthermore, it cannot be applied to samples previously processed using 3' cell barcoding methods since doing so will eliminate the 3' cell barcodes and UMIs, ultimately obscuring single-cell resolution (7). An alternative approach is to couple short reads for gene expression on one sequencing platform (e.g., Illumina) with long reads from another (e.g., Nanopore) to measure phenotype and genotype respectively, as in RAGE-seq, but this approach requires two sequencing platforms and is difficult to scale (7).

Here, we report a method for simultaneously sequencing the transcriptome and *TCRα* and *TCRβ* sequences of T cells using Seq-Well S³, a portable, low-cost platform for massively parallel 3' scRNA-seq, in combination with a TCR target enrichment strategy (18). Seq-Well S³ enables confinement of single cells—whether from human, mouse, non-human primates (NHPs), or the like—with barcoded poly(dT) mRNA capture beads in a PDMS array of sub-nanoliter wells sealed with a semi-permeable membrane; this allows efficient cell lysis through fast solution exchange and improves transcript capture (18, 19). In Seq-Well S³, we obtain a 3' barcoded WTA product as outlined in the online Seq-Well S³ protocol (<http://shaleklab.com/wp-content/uploads/2019/07/SeqWell-S3-Protocol.pdf>). The method as described below also works for enriching TCR repertoires from Drop-seq and Slide-Seq libraries, and can be adapted with minor modification to other 3' barcoded WTA products, such as those generated using 10x 3' kits, or, with additional changes, to 5' barcoded WTA products (15, 16, 18).

From the amplified WTA product, a fraction of the material is used to generate a 3' barcoded whole transcriptome scRNA-seq library. Subsequent sequencing enables mRNA expression quantification from thousands of cells at once. Another fraction of the WTA product is used to capture and enrich *TCRα* and *TCRβ* transcripts contained within the WTA product, and then targeted sequencing is performed on those specific transcripts. Collectively, this approach recovers both the full transcriptome and the TCR sequence of each T cell, enabling direct linkage of T cell clonotype to T cell phenotype. While this chapter will focus on TCR transcript enrichment and sequencing, other sequences of interest, such as oncogenes or viral nucleic acids, can be similarly enriched by appropriately adapting the various primers provided (20, 21). We

also note that this method can be applied to multiple species, as evidenced in the notes regarding primer design (*see* **Notes 1, 2**).

More specifically, the single-cell TCR sequencing protocol presented here utilizes biotinylated probes targeting the *TCRα* and *TCRβ* constant regions (Fig. 2). After affinity capture, the remaining WTA product, which is enriched for *TCRα* and *TCRβ* complementary DNAs (cDNAs), is subsequently amplified using the same universal primer site (UPS) employed to obtain the original WTA library. Through primer extension, a second 5' universal primer site (UPS2) linked to the *TCRα* and *TCRβ* variable regions is then used to modify the pool of TCR nucleic acids. Finally, PCR amplification using a pair of primers that bind to the UPS2 site as well as an extended version of the UPS yields an enriched 3' barcoded single-cell TCR sequencing library. This CDR3-enriched library can be sequenced on an Illumina MiSeq (or NextSeq) using custom sequencing primers specific to the constant regions of the *TCR* genes to recover the variable region of the *TCRα* and *TCRβ* chains. Afterward, these single-cell *TCRα* and *TCRβ* sequences can be matched to whole transcriptome profiles obtained via traditional 3' scRNA-seq methods based on shared 3' cell barcodes (7). In short, sequencing reads are aligned to known TCR sequences and the V region, J region, and CDR3 sequences are identified for each UMI using a computational pipeline implemented in Python. This alignment pipeline returns a table containing these annotations and their corresponding bead barcodes. TCR variable regions may then be matched to corresponding cells from the Seq-Well S³ output by matching the TCR output barcodes to those in the corresponding Seq-Well S³ library.

Below, we provide a detailed walk through of both the experimental and computational pipeline for capturing and annotating single-cell *TCRα* and *TCRβ* sequences.

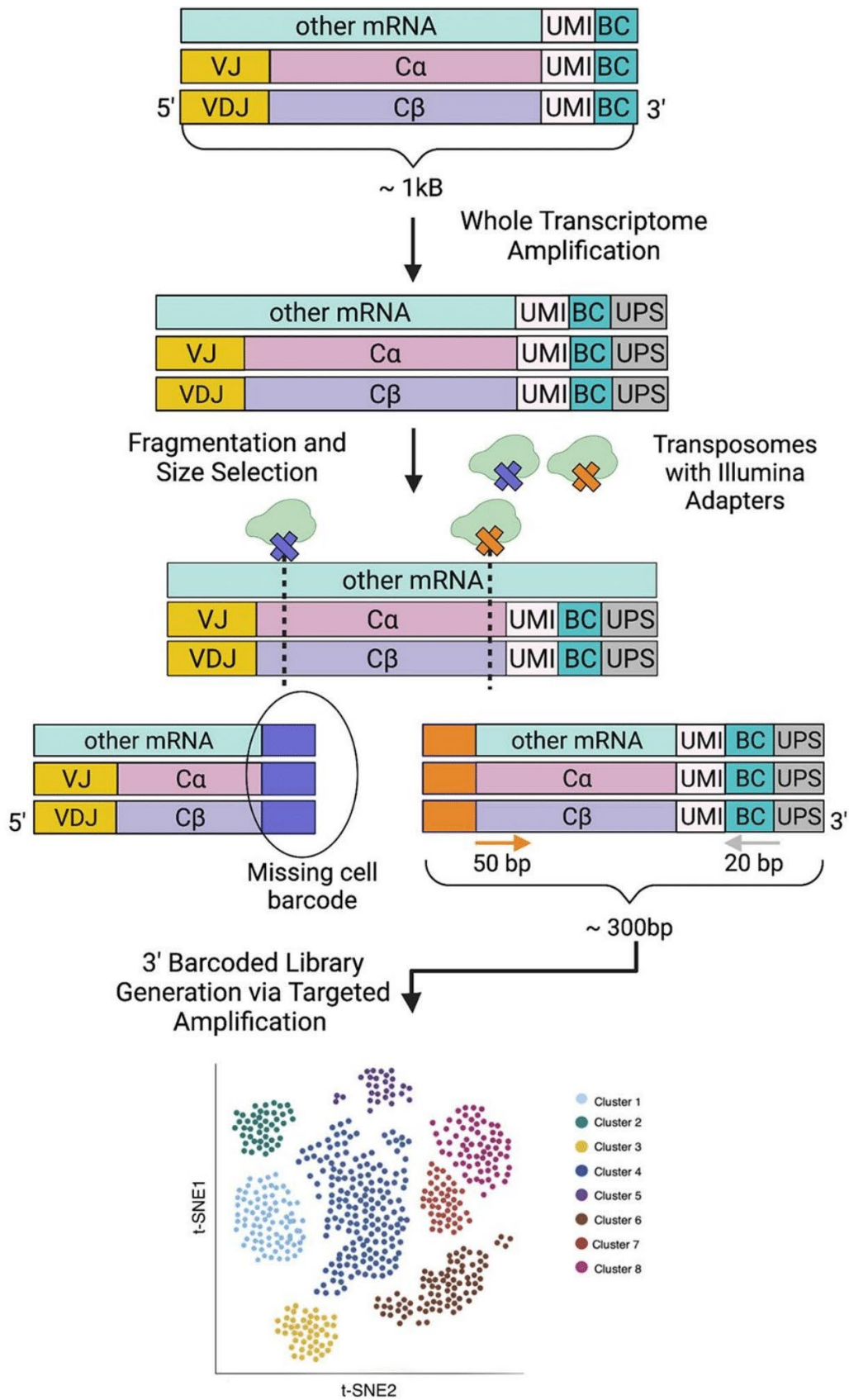


Fig. 1 Using Seq-Well S³, a 3' scRNA-seq method, captured mRNA transcripts are amplified via whole transcriptome amplification (WTA). WTA products, including TCR transcripts, undergo fragmentation, tagmentation, and size selection, resulting in a 3' barcoded library of roughly 300-500 bp size fragments that can be sequenced. For captured TCR transcripts, this process significantly minimizes the chances of obtaining the CDR3 reads in the final library, as this region is approximately 1 kB away from the 3' barcoded region.

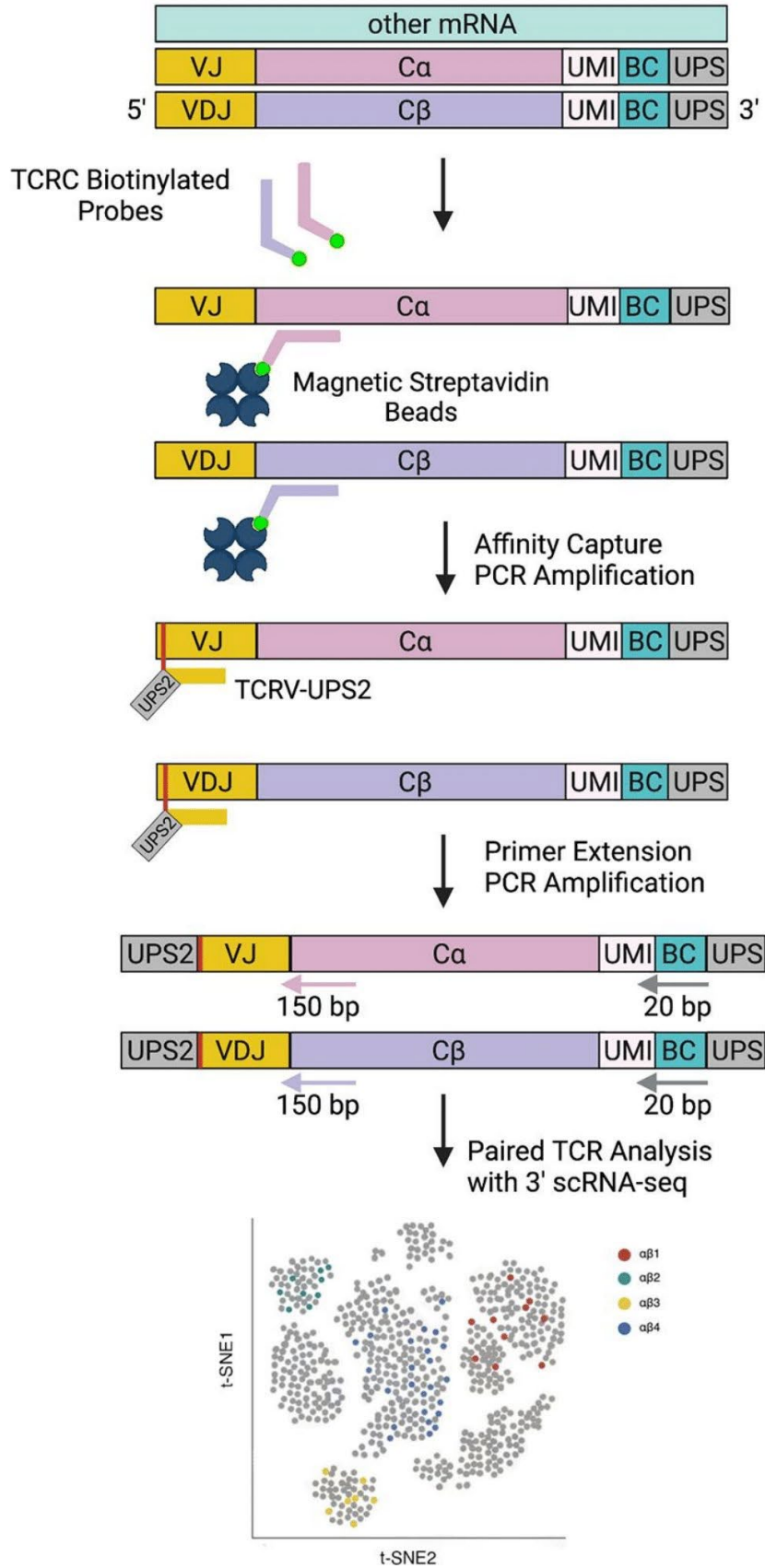


Fig. 2 The amplified WTA product which contains the *TCR α* and *TCR β* transcripts undergoes a 3' cell barcoded library preparation as outlined in the Seq-Well S³ protocol. Affinity capture with species-specific biotinylated oligonucleotides (biotin shown in green) targeting TCR constant regions enriches 3' barcoded TCR transcripts. A sequencing library of these products is then prepared via primer extension with species-specific V-region primers containing UPS2 handles and subsequent PCR amplification with UPS2 and an extended version of UPS. Finally, species-specific custom sequencing primers are used for sequencing the CDR3 region on an Illumina MiSeq and analyzed together with the full transcriptomic data.

2. Materials

2.1. Seq-Well S³ TCR Recovery Materials

1. xGen Hybridization and Wash Kit (IDT, USA), store at 4 °C and -20 °C respectively.
2. Biotinylated oligonucleotide probes for TCR enrichment, store at -20 °C (*see Note 2*).
3. UPS, stock concentration is 100 μ M, store at -20 °C.

5'AAGCAGTGGTATCAACGCAGAGT

4. Nuclease-free water.
5. DNA Lo-Bind tubes, 1.5 mL.
6. PCR strip tubes.
7. 15 mL conical tube.

8. Pipettes (P2, P20, P200, P1000).
9. Low retention, sterile, filtered pipet tips (20 μ L, 200 μ L, 1000 μ L).
10. Kapa Hifi Hotstart Readymix (2 \times), store at -20°C.
11. AMPure XP beads (Beckman Coulter, USA), store at 4 °C.
12. Reagent reservoirs.
13. Ethanol, pure (200 proof, anhydrous).
14. High Sensitivity D5000 Reagents (Agilent, USA), store at 4 °C.
15. High Sensitivity D5000 ScreenTape (Agilent, USA), store at 4 °C.
16. TCRV-UPS2 TCR α (*see Note 1* for primer design).
17. TCRV-UPS2 TCR β .
18. 96-well plate.
19. UPS-mod-N50x, stock concentration is 1 mM, store at -20 °C (*see Note 3* for option for using Illumina barcodes).

5'AATGATACGGCGACCACCGAGATCTACACGCCTGTCCGCGGA
AGCAGTGGTATCAACGCAGAGT*A*C

20. UPS-N70x, stock concentration is 1 mM, store at -20 °C.

5'CAAGCAGAAGACGGCATAACGAGATGTCTCGTGGGCTCGG

21. Custom sequencing primer, stock concentration is 100 μ M, store at -20 °C (*see Note 4*).

22. Seq-Well S³ Sequencing Primer, stock concentration is 100 μM, store at -20 °C.

5'GCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGTAC

23. Illumina MiSeq 150 cycle kit.

2.2. Equipment

1. Thermocycler.
2. Vortex.
3. Magnet for PCR strip tubes.
4. DynaMag-2 Magnet.
5. Thermomixer.
6. Mini centrifuge.
7. Centrifuge.
8. Plate/PCR strip vortex.
9. Qubit Fluorometer.
10. Kapa Library Quantification qPCR.
11. Advanced Analytical Fragment Analyzer.
12. Illumina MiSeq.

2.3. Computer

A computer capable of running Docker and access to Google Cloud.

2.4. Data Matrix Files

Digital expression matrices from Seq-Well S³ libraries including cell barcodes for each library.

3. Methods

Carry out all temperature sensitive procedures at suitable temperature as illustrated (*see Note 5*).

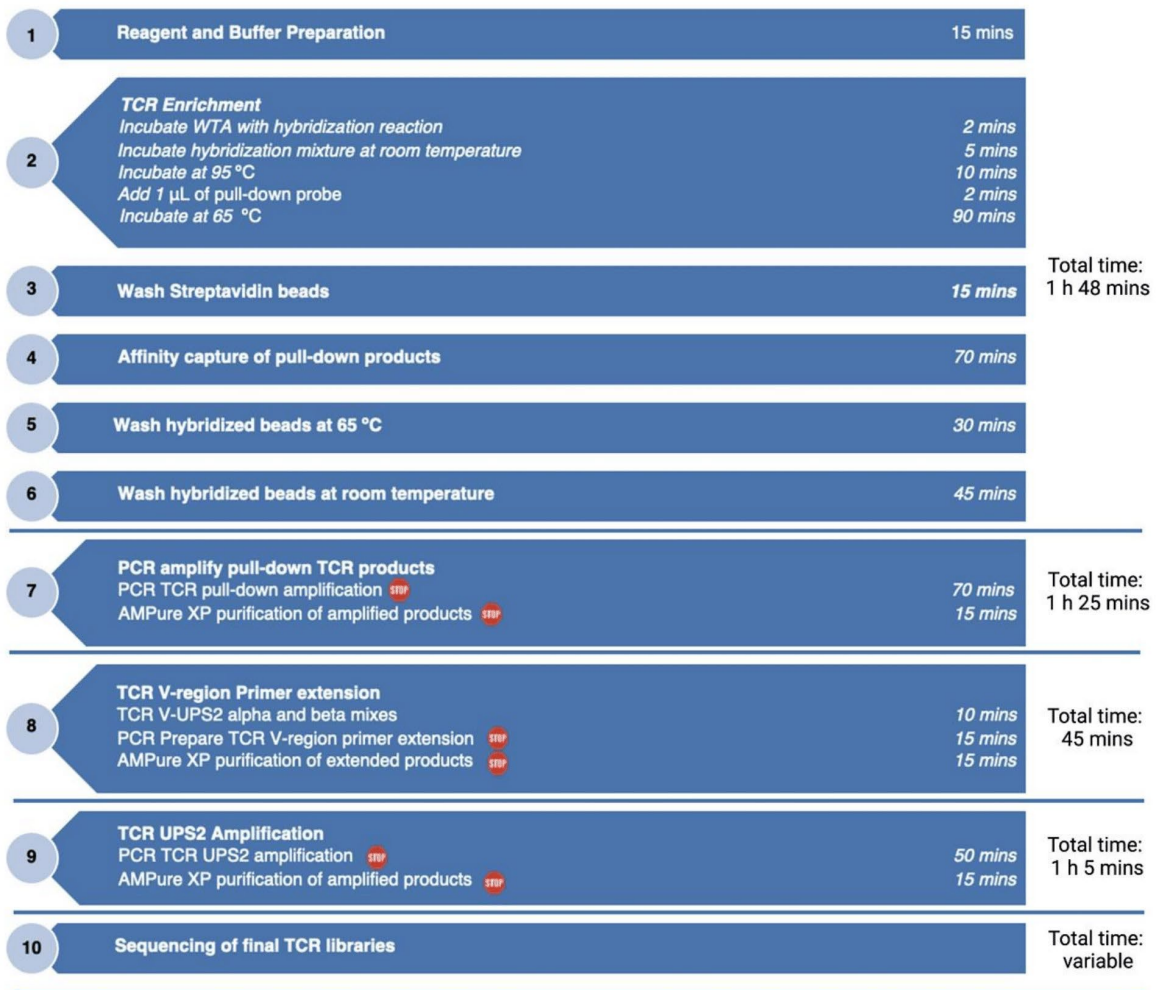


Fig. 3 Workflow depicting the full TCR enrichment protocol with total time of each step specified as well as stopping points in the protocol.

3.1. TCR Enrichment from 3' Barcoded Single-cell Libraries

3.1.1. Reagent Preparation

1. Dilute Wash Buffer 1, Wash Buffer 2, Wash Buffer 3, and Stringent Wash Buffer to working concentrations. Aliquot in 0.2 mL PCR strip tubes and store at 4 °C.
2. Dilute Bead Washing Buffer to working concentration. Aliquot in a 50 mL conical tube and store at 4 °C.
3. Allow Wash Buffer 1, Wash Buffer 2, Wash Buffer 3, Stringent Wash Buffer, and Bead Washing Buffer to come to room temperature before beginning (*see Note 5*).

3.1.2. TCR Enrichment

1. TCR enrichment is performed with whole transcriptome amplified (WTA) products which are generated by performing Seq-Well S³/Drop-seq/Slide-seq methods.
2. Set up a 12.5 µL hybridization reaction with 3.5 µL of WTA product:

Master Mix 1 (per sample/reaction)	1x (μL)
UPS (50 μM)	0.8
Cot-1 DNA	0.5
2× Hybridization Buffer	8.5
xGen Buffer Enhancer	2.7
Total volume	12.5

Table 1 Reagents necessary for master mix 1 of TCR enrichment.

3. Mix 3.5 μL of WTA product with the hybridization reaction mix in 0.2 mL PCR strip tubes.
4. Incubate the tube strip at room temperature for 5 mins. Then, incubate at 95 °C for 10 mins (*see Note 6*).
5. Remove mixture to room temperature and add 1 μL of pull-down probe mix (*see Note 2*). Vortex and centrifuge the strip tubes to mix and collect the liquid respectively.
6. Incubate the mixture at 65 °C for 1 hr to allow for hybridization of pull-down probes to the WTA libraries.

3.1.3. Wash Streptavidin Beads

1. Aliquot 50 μL of streptavidin beads per sample into a 1.5 mL microcentrifuge tube. Place the tube on a DynaMag magnetic stand and allow for the beads to pellet.

2. Remove supernatant and add equal volume of bead wash buffer (BWB) into the tube. (*see Section 3.1.1*)
3. Vortex the tube to mix, and place back on the magnetic stand for beads to pellet.
4. Repeat steps 2-3
5. Take the tube off the magnet and add an equal volume of bead wash buffer.
6. Aliquot 50 μ L (per sample) of the mixture into PCR strip tubes.
7. Place the tubes on a PCR strip magnet and remove the supernatant once the hybridized mixture from **Section 3.1.2** is ready (*see Note 7*).

3.1.4. Affinity Capture of Pull-down Products

1. Add the hybridized mixture from **Section 3.1.2** into the prepared streptavidin beads.
2. Vortex and spin down briefly to mix and collect liquid respectively.
3. Incubate mixture at 65 °C for 45 min. Intermittently vortex the mixture every 10 min to ensure the beads are suspended in the mixture (*see Note 8*).
4. During the incubation, preheat aliquoted Wash Buffer 1 (WB1) and Stringent Wash Buffer (SWB) at 65 °C (*see Note 9*). Preheat at

least 100 μL /sample of WB1 and 400 μL /sample of SWB (2x200 μL).

3.1.5. Wash hybridized beads at 65 degrees

1. While the strip tubes are at 65 °C, add 100 μL of heated WB1 into each tube, pipette up and down gently, scraping the beads from the side of the tubes.
2. Place the PCR strip tube on the magnet, aspirate supernatant, and take tubes off the magnetic stand and return to 65 °C (*see Note 10*).
3. While the strip tubes are at 65 °C, add 200 μL of heated SWB into each tube, and pipette up and down gently resuspend the bead pellets.
4. Incubate mixture at 65 °C for 5 min.
5. Repeat steps 2-4.
6. Aspirate final wash of SWB.

3.1.6. Wash hybridized beads at Room Temperature

1. Add 200 μL of room temperature WB1 into each reaction, pipette up and down to resuspend beads.
2. Vortex for 2 min to mix, spin briefly to collect liquid, and place on a magnetic stand.
3. Aspirate and add 200 μL of room temperature WB2, pipette up and down to resuspend beads.

4. Vortex for 1 min to mix, spin to collect liquid, and place on a magnetic stand.
5. Aspirate and add 200 μL of room temperature WB3, pipette up and down to resuspend beads.
6. Vortex for 30 seconds to mix, spin to collect liquid, and place on a magnetic stand.
7. Aspirate and resuspend the beads in 20 μL of nuclease-free water.

3.1.7. PCR Amplify Pull-down TCR Products

1. Set up a 25 μL of PCR reaction with hybridized mixture. For each mixture, perform 5 PCR reactions, using 2 μL of mixture in each reaction (using a total of 10 out of 20 μL):

TCR pull-down amplification PCR system	1x (μL)
Kapa Hifi Hotstart Readymix (2 \times)	12.5
UPS (10 μM)	2
Nuclease-free water	8.5
Hybridized mixture	2.5
Total volume	25

Table 2 Reagents necessary for TCR pull-down amplification PCR system.

2. Run the PCR using the following conditions:

Cycles	Temperature (°C)	Time
1	95	3 minutes
25	98	40 seconds
	67	20 seconds
	72	1 minute
1	72	5 minutes

Table 3 Outline of PCR program for TCR pull-down amplification.

3. After amplification, pool all 5 PCR reactions (20 μ L for each sample) into a single tube or well in a 96 well plate (for 125 μ L total).

3.1.8. Purify TCR pull-down products and assess quality

1. Add 32 μ L AMPure XP reagent to pooled PCR product from previous step (0.64 \times volumetric ratio); mix by pipetting (*see Note 11*).
2. Allow the tubes to incubate for 5 min at room temperature.
3. Place the tube on the rack for 3 min and remove the supernatant.
4. Add 100 μ L of fresh 80% ethanol to the tube. After addition, shift the tube on the magnet 6 times to allow beads to pass through the ethanol solution to the opposite side of the tube.
5. Remove 80% ethanol wash solution.
6. Repeat steps 3-5 twice.

7. After the third ethanol wash, remove any residual ethanol and allow beads to air dry for 4 min (*see Note 12*).
8. Resuspend beads in 15 μ L nuclease-free water.
9. Place the tubes on the magnet stand and transfer the 15 μ L supernatant, which contains the eluted full-length TCR alpha and beta products, to a fresh 96-well plate.
10. Assess quality of the enrichment via Advanced Analytical Fragment Analyzer or running an Agilent D5000 High Sensitivity Screentape (*see Note 13*) (Fig. 4).
11. Determine concentration using Qubit Fluorometer or Tape report.

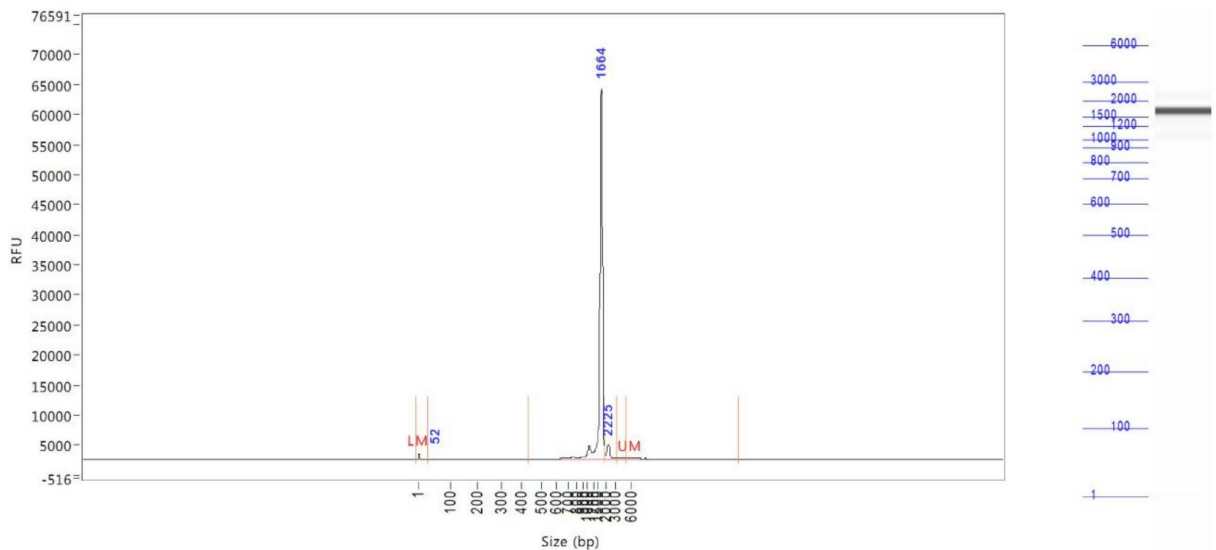


Fig. 4 Expected size distribution of amplified TCR α and TCR β pull-down products (*see Note 14*).

3.2. TCR V-region Primer Extension

1. Prepare separate TCRV-UPS2 alpha and beta mixes by equimolar pool each of the primers and dilute as necessary to 10 μM (total concentration) (*see Note 15*).
2. Set up a 25 μL of PCR reaction with enriched full-length TCR products.
Make two separate reactions (one for TCR α and one for TCR β chain).

TCR V-region primer extension PCR system	1x (μL)
Kapa Hifi Hotstart Readymix (2 \times)	12.5
TCRV-UPS2 TCR α and TCR β (10 μM)	2.5
Nuclease-free water	6
Hybridized mixture	4
Total volume	25

Table 4 Reagents necessary for TCR V-region primer extension PCR system.

3. Perform the primer extension using the following conditions:

Cycles	Temperature ($^{\circ}\text{C}$)	Time
1	95	5 minutes
1	55	30 seconds
1	72	2 minutes

Table 5 Outline of PCR program for TCR V-region primer extension.

4. After extension is complete, add 25 μL of water into each reaction, to bring the total to 50 μL .
5. Purify using AMPure XP beads as previously described in 3.1.8. Elute in 15 μL of water.

3.3. TCR UPS2 amplification:

1. Set up a 25 μL of PCR reaction with 2.5 μL of TCR primer extension product from 3.2. Split each sample into 4 PCR reactions (*see Note 16*).

TCR UPS2 amplification PCR system	1x (μL)
Kapa Hifi Hotstart Readymix (2 \times)	12.5
UPS-mod-N50x (10 μM)	0.5
UPS2-N70x (10 μM)	0.5
Nuclease-free water	9
TCR primer extension product	2.5
Total volume	25

Table 6 Reagents necessary for TCR UPS2 amplification PCR system.

2. Run the PCR using the following conditions (*see Note 17*):

Cycles	Temperature (°C)	Time
1	95	2 minutes
9-18 cycles	95	30 seconds
	60	30 seconds
	72	1.5 minutes
1	72	5 minutes

Table 7 Outline of PCR program for TCR UPS2 amplification.

Concentration of products from 3.1. (ng/ μ L)	Recommended cycle number for TCR β	Recommended cycle number for TCR α
< 1	16-18	17-18
1-10	14	15
10-20	13	14
20-30	12	13
30-40	11	12
>40	9-10	9-11

Table 8 Recommended number of cycles for TCR UPS2 amplification based on concentration estimation of amplified pull-down products. Precise cycles will need to be tested and adjusted for each sample.

3. Take 12.5 μ L of each reaction and pool for each sample for a final volume of 50 μ L per sample.

- Purify using AMPure XP beads as previously described in 3.1.8. Elute in 15 μ L of water (*see Note 18*).
- Assess size distribution via Advanced Analytical Fragment Analyzer or running an Agilent D5000 High Sensitivity Screentape (Fig. 5).
- Determine final concentration using Kapa Library Quantification (*see Note 19*).

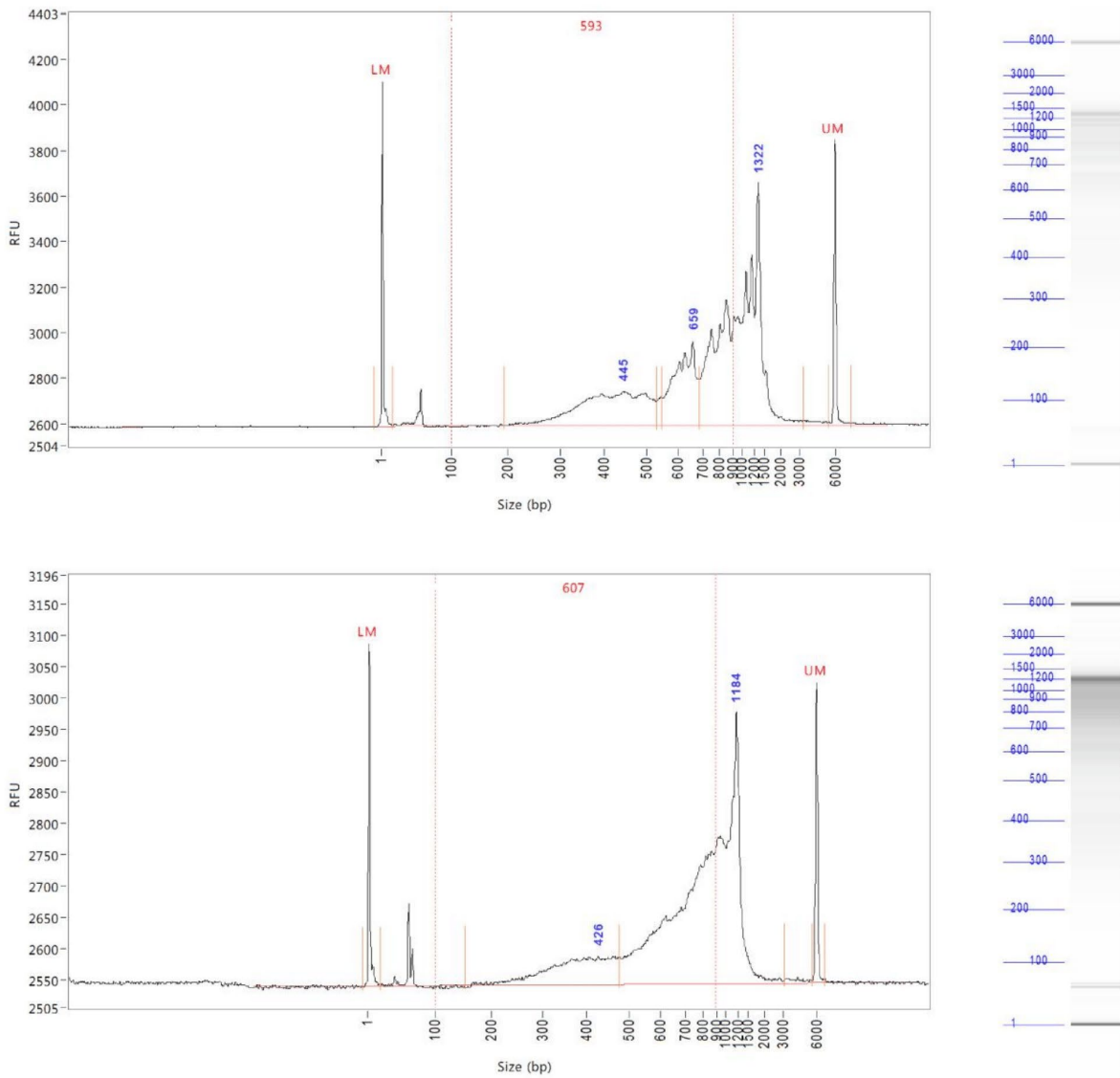


Fig. 5 Expected size distribution of UPS2 amplified TCR α (top) and TCR β products (bottom). TCR α products are expected to be longer than TCR β products (*see Note 14*).

3.4. Sequencing of final TCR libraries:

Final TCR α and TCR β libraries can be pooled based on qPCR concentrations and sequenced on Illumina Miseq using a 150-cycle kit (*see Note 20*). TCR α and TCR β libraries may be pooled or sequenced separately. Sequencing primers are used at a final concentration of 2.5 μ M (*see Note 4*). Species-specific TCR custom sequencing primers are used for Read 1 and the Seq-Well S³ Sequencing Primer is used for index 1. We aim for 8-12 $\times 10^6$ pass filter reads per lane, a cluster density of roughly 450,000/mm² (*see Note 21*). Based on the whole-transcriptome data, we allot ~6,000 T cells per lane.

The sequencing structure is as follows (*see Note 22*):

Read 1: 150 bp

Index 1: 20 bp

3.5. Analysis

After sequencing, computational analysis is required to identify TCRs and connect these with their corresponding cell barcodes. If no Illumina indices are used, the final FASTQ will contain TCR sequences in Read 1 and the cell barcode

with UMI in Index 1. The Illumina bcl2fastq software may be run without barcodes in the sample sheet and a single-ended FASTQ file for the sequenced pool will be output. If Illumina indices are used, see **Section 3.5.1** for demultiplexing instructions. Once FASTQ files are obtained, the TCR alignment and CDR3 identification pipeline should be run (*see Section 3.5.2*).

3.5.1. Demultiplexing with Illumina indices

When Illumina indices are included in the UPS-mod-N50x primer, a sample sheet with these indices should be used as input to bcl2fastq (https://github.com/ShalekLab/tcrgo/blob/master/examples/SampleSheet_indexed_run_example.csv). When running bcl2fastq to generate outputs for the TCRGO pipeline, it is recommended to generate paired-end FASTQ files. This can be done with the ‘--create-fastqs-for-index-reads’ parameter and with the ‘--use-bases-mask’ set to ‘Y140,Y20,I8’ in bcl2fastq. Note that the output of this run will be paired ended fastqs for each barcoded sample where Read 1 contains the 140 bp TCR read and read 2 contains the 20 bp index 2 read with the Seq-Well S³ bead barcode and UMI. At this stage, both FASTQ files represent the reverse complement of the sequences (Fig. 7A), and this format is the expected input of the following pipeline.

3.5.2. TCR alignment and CDR3 identification pipeline

The computational Seq-Well S³ TCR alignment pipeline, TCRGO, identifies and quantifies the TCR regions associated with each cell barcode in the sequencing data resulting from TCR enrichment protocol described in Tu et al. 2019 (7). This pipeline is available at www.github.com/ShalekLab/tcrgo. We have reengineered the Seq-Well S³ TCR recovery pipeline originally presented by Tu et al. to make significant enhancements to the user and developer experience while delivering comparable results. Chiefly, we have improved the performance, reporting, tunability, portability, and maintainability of the computational pipeline in a new Python program, TCRGO. This pipeline consists of four main steps: 1. preprocessing and alignment; 2. filtering using alignment information; 3. CDR3 recovery; and 4. summarizing the data (Fig. 6).

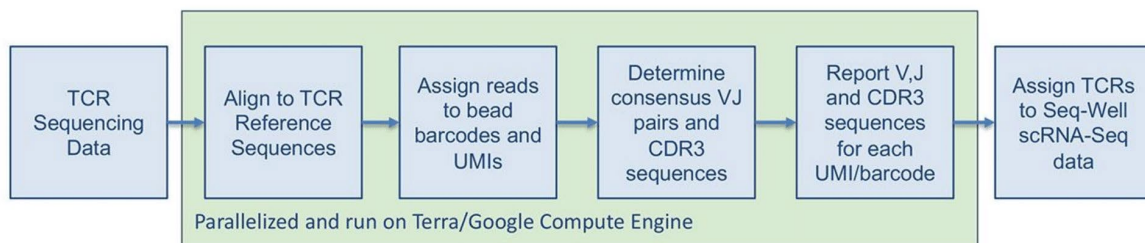


Fig. 6 Overview of TCRGO workflow. Pipeline can be run locally or run on remote servers.

Boxed section is parallelized and optimized to run on the Terra.bio platform.

Running the pipeline

The TCRGO pipeline is available as a Docker image, hosted on Dockstore at <https://dockstore.org/workflows/github.com/ShalekLab/tcrgo/TCRGO>, thus is quickly deployable in cloud computing environments. This pipeline is specifically optimized for use on the Terra.bio platform, a wrapper around Google Compute Engine for computational biology pipelines. To run this pipeline on the Terra platform, FASTQ files should be uploaded to Terra and added to a Terra data table. The pipeline is also available as a series of scripts runnable locally through Python 3.8 with several well-maintained programs as dependencies, including Bowtie 2, Samtools, Picard, and Picard extension Drop-seq Tools (15, 22–24)

Input and reference file formats

The required input files to TCRGO include the TCR sequencing data and a reference FASTA file. Example reference FASTA files for human (*Homo sapiens*) and mouse (*Mus musculus*) are available (<https://github.com/ShalekLab/tcrgo/tree/master/references>). These reference files are multiple sequence FASTA files, which include a separate sequence for each known TCR V and J region. For each region, the description line is written as '>TR[AB][VJC][1-9-]+*' with the following line containing the sequence.

TCR sequencing data may be represented as a single-ended FASTQ file of the 150 bp TCR read and the Seq-Well S³ bead barcode and UMI represented as the index sequence portion of the sequence description (Fig. 7B). Alternatively, the TCR sequencing data may be represented as paired-end FASTQ files with the TCR read on one end and the Seq-Well S³ bead barcode and UMI on the other end (Fig. 7A). When files of this format are input into the TCRGO pipeline, they should be specified with the barcode FASTQ as the first input FASTQ and the TCR FASTQ as the second input FASTQ.

A. tcrseq_R2.fq

```
@M00578:762:JLPRW:1:2107:9088:4212 2:N:0:GAGTAAGA
TAACAACGACAGATCGGCG
+
-A,C@@<:@,C+FFFCFCG@
```

tcrseq_R1.fq

```
@M00578:762:JLPRW:1:2107:9088:4212 1:N:0:GAGTAAGA
ACGGATGAACAATAAGACTGGTTCCTCTTCCAAATGTAAGTACGTAGCCTCCTCCT-
GACCTCGCTGACAGCACAGAGATAAGTGGCTGAGTCTTCAGGCTGGGATCCTTTCTGTCTCTTATACACATCTCC
GAGCCCA
+
C9CCCGCGGGFGGGGE<FGGGGGF<FGGGGGGGFFGGA<@FCGGGGGGGGGGDGEF8F-
GGGGGEGGGGGGGGEGFGGGFFGEEFCGFGGGEFGGFGGDGCFDGFGGFFGGGGGGGGGCGFCG<CFGFGC,C==+@
FE
```

B. tcrseq_no_indices_R1.fq

```
@M00578:762:JLPRW:1:2107:9088:4212 1:N:0:TAACAACGACAGATCGGCG
ACGGATGAACAATAAGACTGGTTCCTCTTCCAAATGTAAGTACGTAGCCTCCTCCT-
GACCTCGCTGACAGCACAGAGATAAGTGGCTGAGTCTTCAGGCTGGGATCCTTTCTGTCTCTTATACACATCTCC
GAGCCCA
+
C9CCCGCGGGFGGGGE<FGGGGGF<FGGGGGGGFFGGA<@FCGGGGGGGGGGDGEF8F-
GGGGGEGGGGGGGGEGFGGGFFGEEFCGFGGGEFGGFGGDGCFDGFGGFFGGGGGGGGGCGFCG<CFGFGC,C==+@
FE
```

Illumina i5 barcode

SeqWell bead barcode

SeqWell UMI

Aligns to V region

Aligns to J region

Aligns to C region

CDR3 region

Fig. 7 Examples of FASTQ file structure for the two types of libraries. **A)** Files output from a run which includes an Illumina barcode after demultiplexing with the ‘create-fastq-for-index-reads’ option. The read 2 FASTQ contains the reverse complement of the UMI and bead barcode while the Read 1 FASTQ contains the reverse complement of the TCR sequence which can be mapped to the V, J, and C regions and the CDR3 can be reconstructed. The Illumina barcode serves to identify the original library of origin of each read in the FASTQ. **B)** Shows the file

output from a run which does not include Illumina barcodes. This file contains the same information as the barcoded files in **A** but relies on the bead barcode to identify the library of origin of the TCR read.

The TCRGO pipeline is highly tunable to user input specifications and allows a variety of input arguments. Several arguments are provided to adjust various thresholds. In the filter queries step, the argument ‘minimum_reads’, with a default of 5, adjusts the minimum threshold of number of reads assigned to a single UMI to include that UMI in the following alignment steps. During the CDR3 recovery step, several thresholds of interest control the outputs. The ‘minimum_frequency’ threshold describes, for a given UMI, the minimum frequency of total reads that must be aligned to the top V/J combination for this UMI to be included. This threshold mostly filters out UMIs with few reads, especially when those reads map to very different places. Relatedly, the ‘minimum_cdr3s’ threshold determines the minimum number of reads mapped to the top V/J combination necessary to proceed with CDR3 assembly for a given UMI. If the library has generally lower read depth and results in lower total reads, another parameter to tune to increase the number of recovered TCRs without relaxing this threshold is the ‘exclude_relatives’ parameter. When selected, this parameter directs the program to use reads mapped to any sequence in the TCR family of the

top V/J regions in assembling the CDR3. For example, if a subset of reads is mapped to TRAV4-1 and another subset to TRAV4-2, reads mapped to both of these regions would be used to construct the CDR3 sequence.

Further input arguments include technical specifications of disc sizes as well as parameters to tune the extent of parallelization of the pipeline. In depth descriptions of these are available at

<https://github.com/ShalekLab/tcrgo#list-of-inputs>.

Preprocessing and alignment

The preprocessing and alignment step formats input files, performs quality control, and aligns the reads to the reference. It prepares the input files for alignment by reformatting them to single-end BAM format to allow the use of the Drop-seq Tools pipeline for quality control. Using components of the Drop-seq Tools program, the reads are tagged with cell and molecular barcodes, low quality barcodes are filtered, and primer and polyA sequences are trimmed. Multimap alignment to V, J, and C TRA and TRB gene segments is performed via Bowtie 2 while all other preprocessing tasks are handled by Samtools, Picard, and Picard extension Drop-seq Tools (15, 22–24). We expect a very high mapping percentage (around 98%) and that a majority of reads map to multiple locations because each read should map to at least a V and J region (Fig. 7A).

Filtering

The filtering step acts to remove reads which do not contain both V and J mapping information. It iterates over the alignment data and distributes reads that multimap to both a V and J segment to lists which can be used in parallel by the CDR3 recovery step.

CDR3 Recovery

The CDR3 recovery step quantifies the alignment identity data for each read, then it resolves the V and J identities for each UMI and determines the UMI's most likely CDR3 sequence using a greedy, optimized, heuristic-driven algorithm.

First, the V and J regions are identified for each read. Often, a single read maps to multiple V or J regions, so the top V and J alignments are determined by Bowtie 2 alignment score, or 'AS' tag. If there is a tie in the alignment score, several metrics are considered. First, the edit distances to the reference sequences are considered. If the tie continues, the Phred quality of the alignment is considered. Then, if the alignments span different lengths of the read, the longer alignment is used. If all options result in a tie, the top alignment is determined alphabetically. If a tie occurs, the information about the tie and the method used to break the tie is recorded for user debugging. Each UMI is assigned the V and J regions to which the most reads were mapped. Users may specify a

minimum frequency threshold to filter for UMIs for which a minimum proportion of reads map to their most frequent V and J regions. This threshold allows for more stringent filtering of UMIs without a conclusive mapping location due to sequencing errors or other technical artifacts.

CDR3 positions file

The `--cdr3-positions-file` argument can be used to input a 2-column TSV file containing the names of FASTA TR(A|B)(V|J) segments and that entry's corresponding CDR3 start or ending codon position in the nucleotide sequence. Each entry name in the CDR3 positions file must exactly match its respective entry in the FASTA. If this file is provided, the CDR3 sequences are identified using the information about the mapped V and J regions and the locations in this file. Since small sequencing errors or annotation differences may shift these locations by one base pair, the amino acid sequence of these regions is reported.

If a CDR3 positions file is not included, a greedy algorithm is used to identify this sequence for each read which maps to the most likely V and J regions. Starting at the V alignment ending position and ending at the J alignment starting position, potential CDR3 sequences are scored based on the absence of internal stop codons, starting codon encoding amino acid cysteine and ending codon encoding the amino acid phenylalanine as these are the conserved starting and stopping codons of the CDR3 region.

To resolve the CDR3 sequence from multiple reads with slightly different sequences, the most frequent CDR3 sequence is used. If there is a tie for the most frequent sequence, a consensus sequence is greedily built from the tied sequences. Ties are found most frequently for UMIs with low numbers of reads.

Output format

The alignment identity and CDR3 information is gathered and tabulated during the summary step as a tab-separated value text file. Each unique UMI and barcode pair is output as a single line in the summary table which includes the read counts for that UMI, V and J alignment information, and CDR3 consensus sequence, example available at <https://github.com/ShalekLab/tcrgo/tree/master/examples> (Table S1). If the ‘collapse’ argument is specified, barcode/UMI pairs are collapsed based on similarity of V/J mapping and CDR3 sequence and an additional output file ending with ‘aggrcollapsed_cdr3_info.tsv’ is generated.

3.5.3. Mapping to Seq-Well S³ cells

After reads are assigned to cell barcodes and UMIs, these TCRs must be assigned to T cells found in the full Seq-Well S³ transcriptomic data. To

do this, we obtain a list of cell barcode sequences from the Seq-Well S³ data matrix and identify the most similar by Hamming distance Seq-Well S³ barcode to each TCR barcode. Example python code for this application is available at <https://github.com/ShalekLab/tergo/tree/master/examples>.

4. Notes

1. V region primer design requires access to an annotation of known V region sequences for the species of interest such as those available in IMGT (25). These primer mixes are designed manually such that each V region may be bound by a maximum of one primer to prevent competitive binding. The binding region is preferred to be found near the middle of the V sequence to allow sufficient sequence length to uniquely identify that region during alignment and allow the sequencing read to span the CDR3 region. These sequences are approximately 60 bases long and have a melting point near 70 °C. The first 30 bases of the sequence include the UPS2 sequence, and the remaining bases of the sequence are complementary to a V region. Primer sequences have been optimized for human and mouse (Table S5B, C) and are undergoing testing for crab-eating macaques (Table S5A).
2. This protocol has been optimized for the following TCR probes which are complementary to the alpha and beta constant regions of the corresponding species (Table S2). When multiple probes for a single constant region are listed, they bind to two overlapping regions of the constant region and should be used if sequential enrichment steps are run.

3. There are two options for demultiplexing TCR regions from pooled sequencing libraries. The first simply uses the Seq-Well S³ bead barcodes to match TCRs to their corresponding cells. If profiling a limited number of T cells, it is unlikely that you will have overlapping Seq-Well S³ bead barcodes between two different T cells; thus, in certain instances, you can pool a handful of TCR libraries together into a single sequencing runs without additional indexing barcodes and still map each read to the original Seq-Well S³ array and T cell of origin. Alternatively, each TCR sample can be indexed by adding a barcode into the UPS2-N50x primers (Table S3). If this is done, primers with different index sequences should be used for each sample and no samples tagged with sequences which differ by 2 bp or fewer should be added to the same sequencing lane.
4. This protocol has been optimized for the following custom read 1 sequencing primers (Table S4). These sequencing primers are complementary to the constant region of the TCR alpha or beta and are species-specific unless otherwise noted.
5. Thaw all buffers at RT, except 10X Wash buffer 1 which may need to be heated at 37°C to ensure all salts go into solution and to resuspend particulates.
6. If using a thermocycler, the lid should be set to 105 °C.
7. Ensure beads are resuspended and that they do not dry in the well before the addition of a sample. Small amounts of Bead Wash Buffer will not interfere with downstream binding of the DNA to the beads.
8. Here we recommend using a thermal shaker to facilitate intermittent vortex.
9. The buffers should be heated at 65 °C for at least 15 min before continuing to wash steps.

10. To minimize non-specific binding of off-target DNA sequences to the capture probes, work quickly to prevent the temperature of the hybridized samples from dropping significantly at temperature sensitive steps.
11. 0.6 means the ratio of AMPure XP reagent volume to sample volume. For example, 0.6 of AMPure XP reagent concentration means 60 μ l AMPure XP reagent in 100 μ l sample (100 x 0.6). Ratio of AMPure XP can be adjusted for improved recovery.
12. Beads should appear matte and have a slight cracked appearance. Prevent over-drying of beads as this can lead to significant yield loss.
13. Certain samples may need to be enriched multiple times, especially samples with low gene recovery or low number of α and β T cells. If the amplification was unsuccessful, we recommend decreasing the number of cycles of PCR amplification from 25 to 12 for subsequent rounds of enrichment.
14. Slight variation between samples is acceptable and expected. After affinity purification, we see peaks that range between 1-2 kB. Upon completing the full pulldown, we expect TCR α peaks around \sim 1.3 kB and TCR β peaks \sim 1.1 kB. DNA contamination from miscellaneous products from WTA are expected and may result in traces with multiple peaks. However, these products will not be sequenced.
15. This protocol has been optimized for the following TCR V-region primers (Table S5A-C).

16. Here to minimize barcode swapping, we split into 4 PCR reactions. Only the UPS-500 primers will have a barcode. You may choose to rely on Seq-Well S³ bead barcodes for demultiplexing and forgo the Illumina barcode.
17. The precise number of cycles needs to be tested and adjusted for each sample. **Table 8** is a recommendation based on the concentration estimation of amplified pull-down products from **Section 3.1.8**.
18. If amplification was unsuccessful after AMPure XP purification, the remaining half of the amplified product can be further amplified for 2-4 PCR cycles.
19. It is expected that the concentration of samples obtained from the fragment analyzer will be higher than the qPCR determined concentration due to the presence of miscellaneous DNA products from the original WTA.
20. For libraries consisting of a high number of cells, a NextSeq could be used. For NextSeq, sequencing primers are used at a final concentration of 2.5 μM . We aimed for $150\text{-}200 \times 10^6$ pass filter reads per lane, a cluster density of roughly $100,000/\text{mm}^2$). Based on the whole-transcriptome data, we allotted $\sim 80,000$ T cells per lane.
21. Due to the large fragment length of the final libraries, we aim to load the flow cells at a lower cluster density.
22. For sequencing with Illumina index barcodes, Read 1 should be assigned as 140 bp, Index 1 should be 20 bp, and Index 2 should be assigned as 8 bp.

References

1. Kirsch IR, Watanabe R, O'Malley JT, et al (2015) TCR sequencing facilitates diagnosis and identifies mature T cells as the cell of origin in CTCL. *Sci Transl Med* 7:308ra158
2. Lossius A, Johansen JN, Vartdal F, et al (2014), High-throughput sequencing of TCR repertoires in multiple sclerosis reveals intrathecal enrichment of EBV-reactive CD8 T cells, <http://dx.doi.org/10.1002/eji.201444662>
3. Schrama D, Ritter C, and Becker JC (2017) T cell receptor repertoire usage in cancer as a surrogate marker for immune responses. *Semin Immunopathol* 39:255–268
4. Taussig MJ (1980) Antigen-specific suppressor T cell factors. *Cancer Immunol Immunother* 9:23–30
5. Farmanbar A, Kneller R, and Firouzi S (2019) RNA sequencing identifies clonal structure of T-cell repertoires in patients with adult T-cell leukemia/lymphoma. *NPJ Genom Med* 4:10
6. Rosati E, Dowds CM, Liaskou E, et al (2017) Overview of methodologies for T-cell receptor repertoire analysis. *BMC Biotechnol* 17:61
7. Tu AA, Gierahn TM, Monian B, et al (2019) TCR sequencing paired with massively parallel 3' RNA-seq reveals clonotypic T cell signatures. *Nat Immunol* 20:1692–1699
8. Hwang B, Lee JH, and Bang D (2018) Single-cell RNA sequencing technologies and bioinformatics pipelines. *Exp Mol Med* 50:1–14
9. Sims S, Willberg C, and Klenerman P (2010) MHC-peptide tetramers for the analysis of antigen-specific T cells. *Expert Rev Vaccines* 9:765–774
10. Huang J, Zeng X, Sigal N, et al (2016) Detection, phenotyping, and quantification of antigen-specific T cells using a peptide-MHC dodecamer. *Proc Natl Acad Sci U S A* 113:E1890–7
11. Ranasinghe S, Lamothe PA, Soghoian DZ, et al (2016) Antiviral CD8⁺ T Cells Restricted by Human Leukocyte Antigen Class II Exist during Natural HIV Infection and Exhibit Clonal Expansion. *Immunity* 45:917–930
12. Dash P, McClaren JL, Oguin TH 3rd, et al (2011) Paired analysis of TCR α and TCR β chains at the single-cell level in mice. *J Clin Invest* 121:288–295
13. Tirosh I, Izar B, Prakadan SM, et al (2016) Dissecting the multicellular ecosystem of metastatic melanoma by single-cell RNA-seq. *Science* 352:189–196
14. Kobschull JM and Zador AM (2018) Cellular barcoding: lineage tracing, screening and beyond. *Nat Methods* 15:871–879
15. Macosko EZ, Basu A, Satija R, et al (2015) Highly Parallel Genome-wide Expression Profiling of Individual Cells Using Nanoliter Droplets. *Cell* 161:1202–1214
16. Rodriques SG, Stickels RR, Goeva A, et al (2019) Slide-seq: A scalable technology for measuring genome-wide expression at high spatial resolution. *Science* 363:1463–1467
17. Klein AM, Mazutis L, Akartuna I, et al (2015) Droplet barcoding for single-cell transcriptomics applied to embryonic stem cells. *Cell* 161:1187–1201
18. Gierahn TM, Wadsworth MH II, Hughes TK, et al (2017), Seq-Well: portable, low-cost RNA sequencing of single cells at high throughput, <http://dx.doi.org/10.1038/protex.2017.006a>
19. Hughes TK, Wadsworth MH, Gierahn TM, et al (2019), Highly Efficient, Massively-Parallel Single-Cell RNA-Seq Reveals Cellular States and Molecular Features of Human Skin Pathology, <https://www.biorxiv.org/content/10.1101/689273v1>

20. Galen P van, Hovestadt V, Wadsworth MH II, et al (2019), Single-Cell RNA-Seq Reveals AML Hierarchies Relevant to Disease Progression and Immunity, <http://dx.doi.org/10.1016/j.cell.2019.01.031>
21. Ziegler CGK, Ferreira I, Clayton K, Stuart H, Hunter J, Walker B, Balasz A, Sigal A, Shalek AK (2021), Paired Host/Viral Single-Cell Genomics to Understand Cellular Responses to HIV-1 Infection [unpublished manuscript]
22. Li H, Handsaker B, Wysoker A, et al (2009) The Sequence Alignment/Map format and SAMtools. *Bioinformatics* 25:2078–2079
23. Langmead B and Salzberg SL (2012) Fast gapped-read alignment with Bowtie 2. *Nat Methods* 9:357–359
24. Broad Institute Picard toolkit, <http://broadinstitute.github.io/picard/>
25. Lefranc, M.-P, Lefranc, G (2003) *The T cell receptor FactsBook*, Academic Press, London, UK

Supplementary Table Captions

Supplementary Table 1 Example output of TCRGO pipeline.

Supplementary Table 2 Biotinylated oligonucleotide probes for TCR enrichment for human (*Homo sapiens*), mouse (*Mus musculus*), and crab-eating macaque (*Macaca fascicularis*, abbreviated MacFas), respectively.

Supplementary Table 3 UPS2_N50x primers with Illumina barcodes.

Supplementary Table 4 Optimized custom TCR sequencing primers for human, mouse and crab-eating macaque. Human/macaque primers are universal for human and macaque species.

Supplementary Tables 5A-5C Species specific variable region primers (A) human (B) mouse, (C) crab-eating macaque (MacFas).

TRAV_nReads	TRAJ_nReads	TRAC_nReads	TRBV_nReads	TRBJ_nReads	TRBC_nReads	UNKN_nReads
22	22	0	0	0	0	0
16	16	0	0	0	0	0
7	7	0	0	0	0	0
18	18	0	0	0	0	0
23	23	0	0	0	0	0
80	80	0	0	0	0	0
1000	1000	0	0	0	0	0
10	10	0	0	0	0	0
8	8	0	0	0	0	0
29	29	0	0	0	0	0
56	56	0	0	0	0	0
58	58	0	0	0	0	0
9	9	0	0	0	0	0
50	50	0	0	0	0	0
14	14	0	0	0	0	0
45	45	0	0	0	0	0
0	0	0	5	5	0	0
7	7	0	0	0	0	0
15	15	0	0	0	0	0
14	14	0	0	0	0	0
17	17	0	0	0	0	0
13	13	0	0	0	0	0
140	140	0	0	0	0	2
8	8	0	0	0	0	0
0	0	0	8	8	0	0
0	0	0	18	18	0	0
26	26	0	0	0	0	0
8	8	0	0	0	0	0
23	23	0	0	0	0	0
9	9	0	0	0	0	0
97	97	0	0	0	0	0
7	7	0	0	0	0	0
137	137	0	0	0	0	0
31	31	0	0	0	0	0
139	139	0	0	0	0	0
23	23	0	0	0	0	0
7	7	0	0	0	0	0
63	63	0	0	0	0	0
44	44	0	0	0	0	0
47	47	0	0	0	0	0
6	6	0	0	0	0	0
12	12	0	0	0	0	0
6	6	0	0	0	0	0
7	7	0	0	0	0	0
7	7	0	0	0	0	0
5	5	0	0	0	0	0
5	5	0	0	0	0	0
7	7	0	0	0	0	0
74	74	0	0	0	0	0
1000	1000	0	0	0	0	0
8	8	0	0	0	0	0
12	12	0	0	0	0	0
7	7	0	0	0	0	0
10	10	0	0	0	0	0
5	5	0	0	0	0	0
10	10	0	0	0	0	0
6	6	0	0	0	0	0
7	7	0	0	0	0	0
12	12	0	0	0	0	0
335	335	0	0	0	0	0
15	15	0	0	0	0	0
87	87	0	0	0	0	0
334	334	0	0	0	0	0
6	6	0	0	0	0	0
31	31	0	0	0	0	0
7	7	0	0	0	0	0
76	76	0	0	0	0	0
5	5	0	0	0	0	0
16	16	0	0	0	0	0
175	175	0	0	0	0	0
16	16	0	0	0	0	0
35	35	0	0	0	0	0
89	89	0	0	0	0	0
7	7	0	0	0	0	0
147	147	0	0	0	0	0
7	7	0	0	0	0	0
40	40	0	0	0	0	0
64	64	0	0	0	0	0
58	58	0	0	0	0	0
13	13	0	0	0	0	0
37	37	0	0	0	0	0
24	24	0	0	0	0	0
86	86	0	0	0	0	0
50	50	0	0	0	0	0
7	7	0	0	0	0	0
146	146	0	0	0	0	0
47	47	0	0	0	0	0
25	25	0	0	0	0	0

Supplementary_Table2_TCRC_Probes

Name	Sequence	Vendor/Service
Human TRBC-1	/5BiosG/GTGTTCACCCACCCRAGGTCGCTGTGTTTGAGCCATCAGAAGCAGAGATCTCCACACCCAAAAGGCCACACTGGTGTGCCTGGCCACAGGC	IDT Ultramer/Standard Desalting
Human TRAC-1	/5BiosG/CTGTCTGCCTATTCACCGATTTTGATTCTCAAACAATGTGTCACAAAGTAAGGATTCTGATGTGTATATCACAGACAAAAGTGTGCTAG	IDT Ultramer/Standard Desalting
Mouse TRBC-1	/5BiosG/AGGATCTGAGAAATGTGACTCCACCCAAGGTCTCCTTTGAGCCATCAAAGCAGAGATTGCAAACAACAAAAGGCTACCCTCGTGT	IDT Ultramer/Standard Desalting
Mouse TRAC-1	/5BiosG/ACATCCAGAACCAGAACCTGCTGTGTACCAGTTAAAGATCCTCGGTCTCAGGACAGCACCCCTCTGCCTGTTACCCGACTTTGACTCCC	IDT Ultramer/Standard Desalting
Macaca fascicularis (MacFas) TRBC-1	/5BiosG/GAGGACCTGAAAAAGGTGTTCCACCCAAGTTGCTGTRTTTGAGCCATCAGAAGCAGAGATCTCCACACCCAAAAGGCCACGCTGGTG	IDT Ultramer/Standard Desalting
MacFas TRBC-2	/5BiosG/GCCATCAGAAGCAGAGATCTCCACACCCAAAAGGCCACGCTGGTGTGCCTGGCCACAGGCTTCTACCCGACCACGTTGGAGCTGAGCTG	IDT Ultramer/Standard Desalting
MacFas TRAC-1	/5BiosG/ATATCCAGAACCCTGACCCTGCCGTGTACCAGCTGAGAGGCTCTAAATCCAATGACACCTCTGTCTGCCTATTACTGATTTGATTCTG	IDT Ultramer/Standard Desalting
MacFas TRAC-2	/5BiosG/TAAATCCAATGACACCTCTGTCTGCCTATTACTGATTTGATTCTGTAATGAATGTGTCACAAAGCAAGGATTCTGACGTGCATATCAC	IDT Ultramer/Standard Desalting

Supplementary_Table3_UPS2_N50x_indexed_primers

Well Position	Name	Sequence
A1	P5-TSO_Hybrid_N501	AATGATACGGCGACCACCGAGATCTACACTAGATCGCGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
B1	P5-TSO_Hybrid_N502	AATGATACGGCGACCACCGAGATCTACACCTCTCTATGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
C1	P5-TSO_Hybrid_N503	AATGATACGGCGACCACCGAGATCTACACTATCCTCTGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
D1	P5-TSO_Hybrid_N504	AATGATACGGCGACCACCGAGATCTACACAGAGTAGAGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
E1	P5-TSO_Hybrid_N505	AATGATACGGCGACCACCGAGATCTACACGTAAGGAGGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
F1	P5-TSO_Hybrid_N506	AATGATACGGCGACCACCGAGATCTACACACTGCATAGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
G1	P5-TSO_Hybrid_N507	AATGATACGGCGACCACCGAGATCTACACAAGGAGTAGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
H1	P5-TSO_Hybrid_N508	AATGATACGGCGACCACCGAGATCTACACCTAAGCCTGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
A2	P5-TSO_Hybrid_N517	AATGATACGGCGACCACCGAGATCTACACGCGTAAGAGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
B2	P5-TSO_Hybrid_N510	AATGATACGGCGACCACCGAGATCTACACCGTCTAATGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
C2	P5-TSO_Hybrid_N511	AATGATACGGCGACCACCGAGATCTACACTCTCTCCGGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
D2	P5-TSO_Hybrid_N513	AATGATACGGCGACCACCGAGATCTACACTCGACTAGGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
E2	P5-TSO_Hybrid_N515	AATGATACGGCGACCACCGAGATCTACACTTCTAGCTGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
F2	P5-TSO_Hybrid_N516	AATGATACGGCGACCACCGAGATCTACACCCTAGAGTGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
G2	P5-TSO_Hybrid_N518	AATGATACGGCGACCACCGAGATCTACACCTATTAAGGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
H2	P5-TSO_Hybrid_N520	AATGATACGGCGACCACCGAGATCTACACAAGGCTATGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
A3	P5-TSO_Hybrid_N522	AATGATACGGCGACCACCGAGATCTACACTTATGCGAGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
B3	P5-TSO_Hybrid_N521	AATGATACGGCGACCACCGAGATCTACACGAGCCTTAGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
C3	P5-TSO_Hybrid_N5_BC01	AATGATACGGCGACCACCGAGATCTACACAACAATGGGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
E3	P5-TSO_Hybrid_N5_BC03	AATGATACGGCGACCACCGAGATCTACACAACCTTGACGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
F3	P5-TSO_Hybrid_N5_BC04	AATGATACGGCGACCACCGAGATCTACACTTCGCTGAGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
H3	P5-TSO_Hybrid_N5_BC06	AATGATACGGCGACCACCGAGATCTACACTTACGCACGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
A4	P5-TSO_Hybrid_N5_BC07	AATGATACGGCGACCACCGAGATCTACACAATGTTCTGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
B4	P5-TSO_Hybrid_N5_BC08	AATGATACGGCGACCACCGAGATCTACACTGTCGGATGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
C4	P5-TSO_Hybrid_N5_BC09	AATGATACGGCGACCACCGAGATCTACACACAGGTATGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
D4	P5-TSO_Hybrid_N5_BC10	AATGATACGGCGACCACCGAGATCTACACTGTAACCTGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
E4	P5-TSO_Hybrid_N5_BC11	AATGATACGGCGACCACCGAGATCTACACACCAACTGGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
F4	P5-TSO_Hybrid_N5_BC12	AATGATACGGCGACCACCGAGATCTACACTGCTCGACGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
H4	P5-TSO_Hybrid_N5_BC14	AATGATACGGCGACCACCGAGATCTACACTCTGGCGAGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
A5	P5-TSO_Hybrid_N5_BC15	AATGATACGGCGACCACCGAGATCTACACAGCATGGAGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
B5	P5-TSO_Hybrid_N5_BC16	AATGATACGGCGACCACCGAGATCTACACTCTCGGTGCGCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
C5	P5-TSO_Hybrid_N5_BC17	AATGATACGGCGACCACCGAGATCTACACAGGTAAGGGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
D5	P5-TSO_Hybrid_N5_BC18	AATGATACGGCGACCACCGAGATCTACACTCGCTAGAGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
E5	P5-TSO_Hybrid_N5_BC19	AATGATACGGCGACCACCGAGATCTACACAGGTGCGAGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
F5	P5-TSO_Hybrid_N5_BC20	AATGATACGGCGACCACCGAGATCTACACTCCTTGGTGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
G5	P5-TSO_Hybrid_N5_BC21	AATGATACGGCGACCACCGAGATCTACACAGTTGCTTGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
H5	P5-TSO_Hybrid_N5_BC22	AATGATACGGCGACCACCGAGATCTACACTATCTGCCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
A6	P5-TSO_Hybrid_N5_BC23	AATGATACGGCGACCACCGAGATCTACACATTATCAAGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
B6	P5-TSO_Hybrid_N5_BC24	AATGATACGGCGACCACCGAGATCTACACTACTTAGCGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
D6	P5-TSO_Hybrid_N5_BC26	AATGATACGGCGACCACCGAGATCTACACTAAGCACAGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
E6	P5-TSO_Hybrid_N5_BC27	AATGATACGGCGACCACCGAGATCTACACATTGTCTGGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
F6	P5-TSO_Hybrid_N5_BC28	AATGATACGGCGACCACCGAGATCTACACGTCGCCTTGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
G6	P5-TSO_Hybrid_N5_BC29	AATGATACGGCGACCACCGAGATCTACACCAATAGTCGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C

H6	P5-TSO_Hybrid_N5_BC30	AATGATACGGCGACCACCGAGATCTACACGTCATCTAGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
A7	P5-TSO_Hybrid_N5_BC31	AATGATACGGCGACCACCGAGATCTACACCAGCAAGGGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
B7	P5-TSO_Hybrid_N5_BC32	AATGATACGGCGACCACCGAGATCTACACGTAACATCGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
C7	P5-TSO_Hybrid_N5_BC33	AATGATACGGCGACCACCGAGATCTACACCAGGAGCCGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
D7	P5-TSO_Hybrid_N5_BC34	AATGATACGGCGACCACCGAGATCTACACGCCTAGCCGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
E7	P5-TSO_Hybrid_N5_BC35	AATGATACGGCGACCACCGAGATCTACACCATGATCGGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
F7	P5-TSO_Hybrid_N5_BC36	AATGATACGGCGACCACCGAGATCTACACGCCGCAACGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
G7	P5-TSO_Hybrid_N5_BC37	AATGATACGGCGACCACCGAGATCTACACCCAACATTGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C
H7	P5-TSO_Hybrid_N5_BC38	AATGATACGGCGACCACCGAGATCTACACGCACATCTGCCTGTCCGCGGAAGCAGTGGTATCAACGCAGAGT*A*C

Supplementary_Table4_TCR_SequencingPrimers

Name	Sequence	Vendor/Service
Human Alpha Sequencing Primer	AGAGTCTCTCAGCTGGTACACGGCAGGGTCAGGITTCTGGATAT	IDT/HPLC
Human Beta Sequencing Primer	CAAACACAGCGACCTCGGGTGGGAACACSTTKTTCAGGTCCT	IDT/HPLC
Mouse Alpha Sequencing Primer	GTCCTGAGACCGAGGATCTTTTAACTGGTACACAGCAGGTTCTGGGTTCTGGATGT	IDT/HPLC
Mouse Beta Sequencing Primer	TGCTTTTGATGGCTCAAACAAGGAGACCTTGGGTGGAGTCACATTTCTCAGATCCT	IDT/HPLC
Human/Macaque Alpha Sequencing Primer	AGAGICTCTCAGCTGGTACACGGCAGGGTCAGGITTCTGGATAT	IDT/HPLC
Human/Macaque Beta Sequencing Primer	CAAACACAGCIACCTIGGGTGGGAACACSTTKTTCAGGTCCT	IDT/HPLC

Supplementary_Table5A_Human_TCRV_Primers

Sequence Name	Sequence
hTRAV1.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG GTC GTT TTT CTT CAT TCC TTA GTC
hTRAV2.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAC GAT ACA ACA TGA CCT ATG AAC GG
hTRAV3-1.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCT TTG AAG CTG AAT TTA ACA AGA GCC
hTRAV4-1.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCT CCC TGT TTA TCC CTG CCG AC
hTRAV5-1.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAA ACA AGA CCA AAG ACT CAC TGT TC
hTRAV6.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAA GAC TGA AGG TCA CCT TTG ATA CC
hTRAV7.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAC TAA ATG CTA CAT TAC TGA AGA ATG G
hTRAV8.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGC ATC AAC GGT TTT GAG GCT GAA TTT AA
hTRAV9.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGA AAC CAC TTC TTT CCA CTT GGA GAA
hTRAV10.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GTA CAG CAA CTC TGG ATG CAG ACA C
hTRAV12.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGA AGA TGG AAG GTT TAC AGC ACA
hTRAV13-1.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGA CAT TCG TTC AAA TGT GGG CGA A
hTRAV13-2.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGG CAA GGC CAA AGA GTC ACC GT
hTRAV14.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GTC CAG AAG GCA AGA AAA TCC GCC A
hTRAV16.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGC TGA CCT TAA CAA AGG CGA GAC A
hTRAV17.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GTT AAG AGT CAC GCT TGA CAC TTC CA
hTRAV18.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGC AGA GGT TTT CAG GCC AGT CCT
hTRAV19.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GTC CAC CAG TTC CTT CAA CTT CAC C
hTRAV20.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGC CAC ATT AAC AAA GAA GGA AAG CT
hTRAV21.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGC CTC GCT GGA TAA ATC ATC AGG A
hTRAV22.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAC GAC TGT CGC TAC GGA ACG CTA
hTRAV23.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCA CAA TCT CCT TCA ATA AAA GTG CCA
hTRAV24.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAC GAA TAA GTG CCA CTC TTA ATA CCA
hTRAV25.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGT TTG GAG AAG CAA AAA AGA ACA GCT
hTRAV26-1.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCA GAA GAC AGA AAG TCC AGC ACC T
hTRAV26-2.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAT CGC TGA AGA CAG AAA GTC CAG T
hTRAV27.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAC TAA CCT TTC AGT TTG GTG ATG CAA
hTRAV29.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCT TAA ACA AAA GTG CCA AGC ACC TC
hTRAV30.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAA TAT CTG CTT CAT TTA ATG AAA AAA AGC
hTRAV34.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCC AAG TTG GAT GAG AAA AAG CAG CA
hTRAV35.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCT CAG TTT GGT ATA ACC AGA AAG GA
hTRAV36.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGG AAG ACT AAG TAG CAT ATT AGA TAA G
hTRAV38.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCT GTG AAC TTC CAG AAA GCA GCC A
hTRAV39.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCC TCA CTT GAT ACC AAA GCC CGT
hTRAV40.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG GCG GAA ATA TTA AAG ACA AAA ACT C
hTRAV41.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGA TTA ATT GCC ACA ATA AAC ATA CAG G

hTRBV2.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGC CTG ATG GAT CAA ATT TCA CTC TG
hTRBV3-1.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GTC TCA CCT AAA TCT CCA GAC AAA GCT
hTRBV4.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCC TGA ATG CCC CAA CAG CTC TC
hTRBVS-48.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCT CTG AGC TGA ATG TGA ACG CCT
hTRBVS-1.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCG ATT CTC AGG GCG CCA GTT CTC T
hTRBV6-1.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GTG GCT ACA ATG TCT CCA GAT TAA ACA A
hTRBV6-23.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCC CTG ATG GCT ACA ATG TCT CCA GA
hTRBV6-4.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGT GTC TCC AGA GCA AAC ACA GAT GAT T
hTRBV6-56.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGT CTC CAG ATC AAC CAC AGA GGA T
hTRBV6-8.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGT CTC TAG ATT AAA CAC AGA GGA TTT C
hTRBV6-9.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGG CTA CAA TGT ATC CAG ATC AAA CA
hTRBV7-2.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GTC GCT TCT CTG CAG AGA GGA CTG G
hTRBV7-3.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCG GTT CTT TGC AGT CAG GCC TGA
hTRBV7-8.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCC AGT GAT CGC TTC TTT GCA GAA A
hTRBV?-46.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GTC TCC ACT CTG AMG ATC CAG CGC A
hTRBV7-7.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGC AGA GAG GCC TGA GGG ATC CAT
hTRBV7-9.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCT GCA GAG AGG CCT AAG GGA TCT
hTRBV9.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCT CCG CAC AAC AGT TCC CTG ACT T
hTRBV10-13.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCA GAT GGC TAY AGT GTC TCT AGA TCA AA
hTRBV10-2.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGT TGT CTC CAG ATC CAA GAC AGA GAA
hTRBV11.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGC AGA GAG GCT CAA AGG AGT AGA CT
hTRBV12-34.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGC TAA GAT GCC TAA TGC ATC ATT CTC
hTRBV12-5.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCT CAG CAG AGA TGC CTG ATG CAA CT
hTRBV13.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GTC TCA GCT CAA CAG TTC AGT GAC TA
hTRBV14.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGC TGA AAG GAC TGG AGG GAC GTA T
hTRBV15.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGA TAA CTT CCA ATC CAG GAG GCC G
hTRBV16.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGC TAA GTG CCT CCC AAA TTC ACC C
hTRBV18.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGG AAC GAT TTT CTG CTG AAT TTC CCA
hTRBV19.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGG TAC AGC GTC TCT CGG GAG AAG A
hTRBV20-1.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGG ACA AGT TTC TCA TCA ACC ATG CAA
hTRBV24-1.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GTG GAT ACA GTG TCT CTC GAC AGG C
hTRBV25-1.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCA ACA GTC TCC AGA ATA AGG ACG GA
hTRBV27-1.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GTA CAA AGT CTC TCG AAA AGA GAA GAG GA
hTRBV28.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGG GGT ACA GTG TCT CTA GAG AGA
hTRBV29.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GGT TTC CCA TCA GCC GCC CAA ACC TA
hTRBV30.UPS2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GCA GAC CCC AGG ACC GGC AGT TCA T

Supplementary_Table5B_Mouse_TCRV_Primer

Sequence Name	Sequence						
TRAV1	ACACTCTTTCCCTACACGACGCTCTCCGATCTAGGGAACCTTTGCTCGGGTCAAC						
TRAV2	ACACTCTTTCCCTACACGACGCTCTCCGATCTGCCTATGAAGGGCAAGAAGTGAAC						
TRAV3	ACACTCTTTCCCTACACGACGCTCTCCGATCTGTGGAGCAKYGCCCTCCTCAC						
TRAV4_A for 4-1_01 4D-2_01 4-2_01 4-2_02	ACACTCTTTCCCTACACGACGCTCTCCGATCTGGTCTGAGATGCAATTTTCTACCA						
TRAV4_B TRAV4x-3	ACACTCTTTCCCTACACGACGCTCTCCGATCTAACTGTGCAGTGGTTCCTACAGAAT						
TRAV4_C TRAV4x-4	ACACTCTTTCCCTACACGACGCTCTCCGATCTAGAATCCAGGGGCAGCCTCATC						
TRAV5 74 Tm	ACACTCTTTCCCTACACGACGCTCTCCGATCTTCCYTTGGTATAAGCAAGARCCTG						
TRAV5_3P 5x-3 all pseudo	ACACTCTTTCCCTACACGACGCTCTCCGATCTAGAGCTTTCTTTCATCCTAAGAGGC						
TRAV6_1 6-1 6-3 6-4	ACACTCTTTCCCTACACGACGCTCTCCGATCTCCRACTCTTTSTGGTATGTCCAATA						
TRAV6_2 6-7	ACACTCTTTCCCTACACGACGCTCTCCGATCTTACGATACACTGCAACTACTCAGC						
TRAV6_3 6-6	ACACTCTTTCCCTACACGACGCTCTCCGATCTACCCTAATCTTTTCTGGTATGTTCG						
TRAV6_4 rest	ACACTCTTTCCCTACACGACGCTCTCCGATCTCCGGAGAAGTCCACAGCTCCT						
TRAV7_1 7-2 7-3 7-5 7-6 without partials	ACACTCTTTCCCTACACGACGCTCTCCGATCTAAKTRCAGCAGAGCCAGAMTC						
TRAV7_2 the rest without partials	ACACTCTTTCCCTACACGACGCTCTCCGATCTTCTGGTGGTACAGACAGCATTCTG						
TRAV8	ACACTCTTTCCCTACACGACGCTCTCCGATCTGTTCAAATGAGMGAGAGAAGCGCA						
TRAV9_1 70C tm	ACACTCTTTCCCTACACGACGCTCTCCGATCTTATGGTGGATCCATTACCTCTCC						
TRAV9_2 all apart from 9-1	ACACTCTTTCCCTACACGACGCTCTCCGATCTGTACCCGGGCGAGGGGCTGC						
TRAV10	ACACTCTTTCCCTACACGACGCTCTCCGATCTTCAACTGCACTTACACAGATACTGC						
TRAV11	ACACTCTTTCCCTACACGACGCTCTCCGATCTGACAACCCTTAAGGTGGTTCAAAC						
TRAV12_1 all but 12-4_all and 12D-1_3	ACACTCTTTCCCTACACGACGCTCTCCGATCTYCACAGACAACAAGAGMCCGAG						
TRAV12_2 both 12-4	ACACTCTTTCCCTACACGACGCTCTCCGATCTGTTATGCAATGAAATCCATAACCC						
TRAV13	ACACTCTTTCCCTACACGACGCTCTCCGATCTCCTTGGTCTGCAGGAGGSSGA						
TRAV14_1 TRAV14-1 and TRAV14-2	ACACTCTTTCCCTACACGACGCTCTCCGATCTCAGCACTTTTTRACTACTTCCCATG						
TRAV14_3 TRAV14-3	ACACTCTTTCCCTACACGACGCTCTCCGATCTGAGAACAGTGCTTTTGACTACTTCC						
TRAV15_not3	ACACTCTTTCCCTACACGACGCTCTCCGATCTTAGTGGAGAGATGGTTTTTCTATTY						
TRAV15_3	ACACTCTTTCCCTACACGACGCTCTCCGATCTGGACAGAAATTAAGTCCAAAC						
TRAV16	ACACTCTTTCCCTACACGACGCTCTCCGATCTGACAATGGACTGTGTATGAAACC						
TRAV17	ACACTCTTTCCCTACACGACGCTCTCCGATCTCCGTGGACCAGCCTGATGCTC						
TRAV18	ACACTCTTTCCCTACACGACGCTCTCCGATCTAGAGTCCCTCGTTTCTGAGTATCC						
TRAV19	ACACTCTTTCCCTACACGACGCTCTCCGATCTACTATTTGCTGGTACAAAAATACC						
TRAV20	ACACTCTTTCCCTACACGACGCTCTCCGATCTTCAACTGCAGTATACGAATGCAGC						
TRAV21	ACACTCTTTCCCTACACGACGCTCTCCGATCTGGAACGAGTACATCTATTGGTACC						
TRAV22	ACACTCTTTCCCTACACGACGCTCTCCGATCTTGATCTGTAGTTACCAAGACAGAGC						
TRAV23	ACACTCTTTCCCTACACGACGCTCTCCGATCTATTCAAGATCCGCTACAGACAGCC						
TRBV1	ACACTCTTTCCCTACACGACGCTCTCCGATCTGCATCTTGAAGAATTCCAGTATCC						
TRBV2	ACACTCTTTCCCTACACGACGCTCTCCGATCTCCACAATGCTATGTATTGGTATAGAC						
TRBV3	ACACTCTTTCCCTACACGACGCTCTCCGATCTAGCAGATGGAGTTTCTGGTTAATTC						
TRBV4	ACACTCTTTCCCTACACGACGCTCTCCGATCTAAGAAATGCTGAAGATTATGTTTAGC						
TRBV5	ACACTCTTTCCCTACACGACGCTCTCCGATCTAACATCTGGGACATAATGCTATGTAC						
TRBV6	ACACTCTTTCCCTACACGACGCTCTCCGATCTTTTGTACATTTTCATACCAACAGAC						
TRBV7	ACACTCTTTCCCTACACGACGCTCTCCGATCTAATAGTCCCTTCTCTGACACC						
TRBV8	ACACTCTTTCCCTACACGACGCTCTCCGATCTCATGCTTTGGTACCGACAAGATCC						
TRBV9	ACACTCTTTCCCTACACGACGCTCTCCGATCTACTGACAGAACTCAAGATGAGGCC						
TRBV10	ACACTCTTTCCCTACACGACGCTCTCCGATCTCAGGGGTTGAAGCTGATCCATTAC						
TRBV11	ACACTCTTTCCCTACACGACGCTCTCCGATCTCCCTTTACTGGTATCGACAGATC						
TRBV12	ACACTCTTTCCCTACACGACGCTCTCCGATCTAGTCCAACAGTTTGATGACTTACT						
TRBV13	ACACTCTTTCCCTACACGACGCTCTCCGATCTYATGTACTGGTATCGGCAGGACAC						
TRBV14	ACACTCTTTCCCTACACGACGCTCTCCGATCTGTGACCCTATTTCTGACATGATAC						
TRBV15	ACACTCTTTCCCTACACGACGCTCTCCGATCTGAAGTGTGAGCCAGTTTCAGGCC						
TRBV16	ACACTCTTTCCCTACACGACGCTCTCCGATCTGTGACAGGGAAGGACAAGAAGC						
TRBV17	ACACTCTTTCCCTACACGACGCTCTCCGATCTGTTCTCAGACTATGAATCATGATACC						
TRBV18	ACACTCTTTCCCTACACGACGCTCTCCGATCTAACACAGTCATGTTATATGGTATAGTC						
TRBV19	ACACTCTTTCCCTACACGACGCTCTCCGATCTTGAACAGATCTTCAAAAAGGCGATC						
TRBV20	ACACTCTTTCCCTACACGACGCTCTCCGATCTGGAGTGTCAAGCTGTGGGTTTTTC						
TRBV21	ACACTCTTTCCCTACACGACGCTCTCCGATCTTACGTTTATTGGTACTACAAGAAACC						
TRBV22	ACACTCTTTCCCTACACGACGCTCTCCGATCTAGCTCATGACAGAACAGACTTGTTTC						
TRBV23	ACACTCTTTCCCTACACGACGCTCTCCGATCTTAAAGTCTGATTTACTTTTCAGAATC						
TRBV24	ACACTCTTTCCCTACACGACGCTCTCCGATCTGTATCAACAGAACCAGAAACAAGAAC						
TRBV25	ACACTCTTTCCCTACACGACGCTCTCCGATCTTTCTTTGTATTACAGAATCAAGAAGC						
TRBV26	ACACTCTTTCCCTACACGACGCTCTCCGATCTGTATTCTGGTATCAACAAAATAAGAAC						
TRBV27	ACACTCTTTCCCTACACGACGCTCTCCGATCTCAAATGTATTGGTATCAGCAAGACC						
TRBV28	ACACTCTTTCCCTACACGACGCTCTCCGATCTCATAGATGAGACTCATGTCTGAATC						

TRBV29	ACACTCTTCCCTACACGACGCTCTCCGATCTAAACAATGTACTGGTATCGACAAGAC						
TRBV30	ACACTCTTCCCTACACGACGCTCTCCGATCTACAGAGCTTGATGTCATGGCAAC						
TRBV31	ACACTCTTCCCTACACGACGCTCTCCGATCTGCCGAGATCAAGGCTGTGGGC						

Supplementary_Table5C_MacFas_TCRV_Primer

Sequence Name	Sequence																		
macfas_TRA_1	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AGT CTG TGA CCC AGC TTG ACA GCC A																		
macfas_TRA_2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA GAA CAG CTC CCT GCA CAT CAC A																		
macfas_TRA_3	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA GGA TCC CAG CCT GAA GAC TCA G																		
macfas_TRA_4	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA TAC TTT CAA GAA CTG CTT GGA A																		
macfas_TRA_5	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA TCC TTC AGT CTC AAG ATC TCA G																		
macfas_TRA_6	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAC TGC ACT TAC AGC AAC AGT GCT T																		
macfas_TRA_7	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAG CCA CAT ACC GTA AAG AAA CCA C																		
macfas_TRA_8	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AGC TCC AGA TGA AAG ACT CTG CCT C																		
macfas_TRA_9	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG GAT GCT GCT GAG TAC TTC TGT GCT G																		
macfas_TRA_10	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG CAG CAG CTG GAG CAG AGT CCT CGG T																		
macfas_TRA_11	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG GAC TCA GCA ACG TAT TTC TGT GCA A																		
macfas_TRA_12	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG TGG ACT CAG CAG TAT ACT TCT GTG C																		
macfas_TRA_13	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG GTT CAA ATA TGG CTA AGA GGC AAG G																		
macfas_TRA_14	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG CAT ATG GAG CAG GGA GAC AAA CAT T																		
macfas_TRA_15	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAT CTC TTC TGG TAT GTC CAG TAC C																		
macfas_TRA_16	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG CTG TAA GCC AGC ATA ACC ACC ATG T																		
macfas_TRA_17	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG TGT AAC TCT CAA CTG CAG TTA TGA A																		
macfas_TRA_18	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG GAA GCG CCC ATC TTC CTG ATG ATA T																		
macfas_TRA_19	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA GCG ACA GCG TCA CAC TGA ACT G																		
macfas_TRA_20	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AGA ATG GCC TCT CTG ACA ATC ACT G																		
macfas_TRA_21	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG CAA ATG GAG CAG TGA AGC AGG AGG G																		
macfas_TRA_22	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AGA CAT TAT ATG GCT TAC ACT GGT A																		
macfas_TRA_23	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG ATC AGG GTT ATT CTA AGT CAA ATG C																		
macfas_TRA_24	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG GGA GAC AAC TTG GTT CTC AAC TGC A																		
macfas_TRA_25	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG CCT CTC TGA TCA TCA CAG AAG ACA G																		
macfas_TRA_26	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AGT CCA GAA AGG AGG GAT TTC AAT T																		
macfas_TRA_27	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG ACT TCT TCA AGC ACA TTT AAC ACC T																		
macfas_TRA_28	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG CTG GCT GCA ACA GCA TCC AGG AGG A																		
macfas_TRA_29	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAT GCT AAA CAT GTC TCC CTG CAT A																		
macfas_TRA_30	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG ACA TAT ATT CTT TCA AAT ACG GAC C																		
macfas_TRA_31	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA CCA TCT GCC CTT GTG AGC GAC T																		
macfas_TRA_32	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG TCA ACA AGG AGA AGA GGA TCC TCA G																		
macfas_TRA_33	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG ACA AGA CAG CCA AAC ATT TCT CTC T																		
macfas_TRA_34	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG ATC ATG TGG TAC CAA CAG TTT CCC A																		
macfas_TRA_35	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG CCC TGA GAG GGC AGC TCT AAC ATT A																		
macfas_TRA_36	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AGC ACA GCT CGA TAG AGC CAG CCA G																		
macfas_TRA_37	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG ACT CTT AAA CAG AGT TTG TTT CAT A																		
macfas_TRA_38	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG TCT CAG CAG GGA CGA TAC AAC ATG A																		
macfas_TRA_39	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG CTG CTG AAG GTC CTA CAT TCC TGA T																		
macfas_TRA_40	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG TCA TGA CTT TCA GTG AGA ACA CAA A																		
macfas_TRA_41	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG CAC ATT TCC TCC TCC CAG ACC ACA G																		
macfas_TRA_42	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG CAG GTA TCA GAC TCA GCC GTG TAC T																		
macfas_TRA_43	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AGG TGA CGC AGA GTC CTG AGG CCC T																		
macfas_TRA_44	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AGA GTC AAT ATT CAG TGA GCT TCC A																		
macfas_TRAV11-3	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAGTTTGAATATCCCGACCTCTCATC																		
macfas_TRAV14-1	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAGTCTGTACAAGTAGCCAGGAGTGG																		
macfas_TRAV46	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG TCC TTA AGC ATG TGC AGC CAA GAG A																		
macfas_TRAV8-6	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAGGAACGACACGCTGAGTACTTCTGTG																		
macfas_TRB_1	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA CAA ACA GAA AGC AAA GAT GGA T																		
macfas_TRB_2	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA CAA AGA TGG ATT GTA CCC CCG A																		
macfas_TRB_3	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA CAA GAG GAC AGC AAG TGA CAC T																		
macfas_TRB_4	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA CAG CCA CTC TGC AAT GCT ATC C																		
macfas_TRB_5	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA CAT GCA CTG AGT CAA CAC AGT C																		
macfas_TRB_6	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA CCC CAA AGC ACC TGA TCA CAG C																		
macfas_TRB_7	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA CGG ATA TTC AGA GAG GAG ACC T																		
macfas_TRB_8	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA CGT GAG CAT CTG GGA CAT GAT T																		
macfas_TRB_9	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA GAA GAG AAA GGA TGT AGC TCT C																		
macfas_TRB_10	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA GAA GAT CAC TCT GGA ATG TTC T																		
macfas_TRB_11	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA GAA TGT TCT GGT ATC GAC AAG A																		
macfas_TRB_12	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA GAG TTC AAG TTC TTG ATT TCC T																		
macfas_TRB_13	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA GAT TTT AAC AAT GAA GCA GAC A																		
macfas_TRB_14	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA GCC AAG CAG GGA TGT CTG TCA A																		
macfas_TRB_15	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA GGA AGG ACA GAA TGT GAC CCT G																		
macfas_TRB_16	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA GGA ATT CAT TCT ACA GTG TTC C																		
macfas_TRB_17	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA GGA CAC AGT CAT GTT TAT TGG T																		
macfas_TRB_18	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA GGA GTT CTC AGG ATA TGT GCC A																		
macfas_TRB_19	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA TCA GTC TCA GCC CCA TGC CTG G																		
macfas_TRB_20	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA TGA GAG GAC AGC AAG CGA CAT T																		
macfas_TRB_21	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA TGG GAC AAG AAG TAA CTA TGA G																		
macfas_TRB_22	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA TGG GAC AAG CAG TGA CTC TCA G																		
macfas_TRB_23	TCG TGG GCT CGG AGA TGT GTA TAA GAG ACA GAG AAA TTC CAG GTC CTG AAG ACA GGA C																		

